

# Vector Paro / Vector Paro Pro

EN



Installation and operating instructions

CE 0297

9000-615-28/02



 **DÜRR  
DENTAL**

1903V004



# Contents



## Important information

<b>1</b>	<b>About this document</b> . . . . .	3
1.1	Warnings and symbols . . . . .	3
1.2	Copyright information . . . . .	4
<b>2</b>	<b>Safety</b> . . . . .	5
2.1	Intended purpose . . . . .	5
2.2	Intended use . . . . .	5
2.3	Improper use . . . . .	5
2.4	General safety information . . . . .	5
2.5	Specialist personnel . . . . .	5
2.6	Electrical safety . . . . .	6
2.7	Only use original parts . . . . .	6
2.8	Transport . . . . .	6
2.9	Disposal . . . . .	7



## Product description

<b>3</b>	<b>Overview</b> . . . . .	8
3.1	Vector Paro / Vector Paro Pro . . . . .	8
3.2	Scope of delivery . . . . .	10
3.3	Accessories . . . . .	10
3.4	Consumables . . . . .	10
3.5	Wear parts and replacement parts . . . . .	11
<b>4</b>	<b>Technical data</b> . . . . .	12
4.1	Type plate . . . . .	14
4.2	Order number and serial number for the handpieces . . . . .	14
4.3	ID number for tool kits . . . . .	14
4.4	Evaluation of conformity . . . . .	15
<b>5</b>	<b>Operation</b> . . . . .	15
5.1	Handpieces . . . . .	15
5.2	Fluid container . . . . .	15
5.3	Instrument change . . . . .	16



## Assembly

<b>6</b>	<b>Requirements</b> . . . . .	17
<b>7</b>	<b>Installation</b> . . . . .	17
7.1	Establishing the electrical connections . . . . .	17
7.2	Connecting the flexible foot switch . . . . .	18
<b>8</b>	<b>Commissioning</b> . . . . .	19
8.1	Function check . . . . .	19



## Usage

<b>9</b>	<b>Components</b> . . . . .	20
9.1	Paro handpiece . . . . .	20
9.2	Scaler handpiece . . . . .	23
9.3	Instruments and tool kits . . . . .	24
9.4	Steri-box . . . . .	27
9.5	Flexible foot switch . . . . .	27
9.6	Service kit . . . . .	28
9.7	Fluid container . . . . .	28
9.8	Vector Fluid Polish . . . . .	28
9.9	Cleaning components . . . . .	29
9.10	Vector toolcard . . . . .	29
<b>10</b>	<b>Operation</b> . . . . .	30
10.1	Display/handling . . . . .	30
10.2	Adjustment options . . . . .	31
10.3	Preparing the device for treatment . . . . .	32
<b>11</b>	<b>Treatment</b> . . . . .	36
11.1	Preparation . . . . .	36
11.2	Treatment using the Paro handpiece . . . . .	36
11.3	Treatment with a Scaler handpiece . . . . .	38
11.4	Use of the Scaler instruments . . . . .	39
11.5	After every treatment . . . . .	40
<b>12</b>	<b>Cleaning</b> . . . . .	40
12.1	Cleaning of the outside surfaces . . . . .	40

12.2	Activating the cleaning process of the device . . . . .	41	15.9	Changing the support seal . . . . .	64
12.3	Cleaning the fluid container . . . . .	43	15.10	Checking the function of the tool kit cover . . . . .	65
12.4	Cleaning the sleeve part and adapter of the handpiece hose . . . . .	43	15.11	Changing the sealing ring of the resonant ring . . . . .	66
<b>13</b>	<b>Reprocessing . . . . .</b>	<b>44</b>	15.12	Changing the resonant ring . . . . .	67
13.1	Risk analysis and categorisation . . . . .	44	15.13	Changing the interchangeable bushing with union nut . . . . .	67
13.2	Reprocessing procedure in accordance with EN ISO 17664 . . . . .	44			
13.3	Preparation at the operating location . . . . .	46			
13.4	Dismantling the handpiece . . . . .	47			
13.5	Manual cleaning, intermediate rinsing, disinfection, final rinsing, drying in the cleaning and disinfection bath . . . . .	47			
13.6	Manual cleaning, intermediate rinsing, disinfection, final rinse, drying in ultrasonic bath . . . . .	50			
13.7	Automatic cleaning, intermediate rinsing, disinfection, final rinse, drying . . . . .	53			
13.8	Check for function . . . . .	53			
13.9	Packing . . . . .	54			
13.10	Steam sterilising . . . . .	54			
13.11	Issue clearance for the parts for sterilisation . . . . .	54			
13.12	Storing parts for sterilisation . . . . .	55			
<b>14</b>	<b>Treatment breaks for more than 24 hours . . . . .</b>	<b>56</b>			
14.1	Cleaning and disinfecting the hose system . . . . .	56			
14.2	Initial start-up after a break in treatment for more than 24 hours . . . . .	57			
<b>15</b>	<b>Maintenance . . . . .</b>	<b>58</b>			
15.1	Service kit . . . . .	58			
15.2	Changing the seal . . . . .	58			
15.3	Changing the valve in the fluid container . . . . .	59			
15.4	Checking instrument wear . . . . .	59			
15.5	Replacing the light conductor in the Scaler handpiece . . . . .	60			
15.6	Replacing the sealing screw . . . . .	60			
15.7	Inserting or changing the battery in the flexible foot switch . . . . .	61			
15.8	Changing the nozzle in the Paro handpiece . . . . .	63			



**Troubleshooting**

<b>16</b>	<b>Tips for operators and service technicians . . . . .</b>	<b>70</b>
-----------	---	-----------

 **Important information**

## 1 About this document

These installation and operating instructions represent part of the unit.



If the instructions and information in these installation and operating instructions are not followed, Dürr Dental will not be able to offer any warranty or assume any liability for the safe operation and the safe functioning of the unit.

### 1.1 Warnings and symbols

#### Warnings

The warnings in this document are intended to draw your attention to possible injury to persons or damage to machinery.

The following warning symbols are used:



General warning symbol

The warnings are structured as follows:



#### SIGNAL WORD

##### Description of the type and source of danger

Here you will find the possible consequences of ignoring the warning

- Follow these measures to avoid the danger.

The signal word differentiates between four levels of danger:

- **DANGER**  
Immediate danger of severe injury or death
- **WARNING**  
Possible danger of severe injury or death
- **CAUTION**  
Risk of minor injuries
- **NOTICE**  
Risk of extensive material/property damage

#### Other symbols

These symbols are used in the document and on or in the unit:



Note, e.g. specific instructions regarding efficient and cost-effective use of the unit.



Comply with the Operating Instructions.



Protection class II



CE labelling with the number of the notified body



Manufacturer



Health Industry Bar Code (HIBC)



Serial number



Order number



Lot designation



Dispose of correctly in accordance with EU Directive 2012/19/EU (WEEE).



Type BF application part



Steam sterilise at 134 °C



Steam sterilise at 135°C



Not sterile



Recycling



The device contains a battery.



Use suitable tools.



The seal must only be removed by a qualified expert.





On/off switch





Amplitude reduction

 Amplitude increase

 Switch off and de-energise the unit (e.g. unplug from mains).


 Wear protective gloves.


 Wear protective goggles.

 Use a face mask.

 Use protective clothing.

 Rinse with water.

 Rinse with instrument cleaner.

 Rinse with instrument disinfectant.

## 1.2 Copyright information

All circuits, processes, names, software programs and units mentioned in this document are protected by copyright.

The Installation and Operating Instructions must not be copied or reprinted, neither in full nor in part, without written authorisation from Dürer Dental.

## 2 Safety

Dürr Dental has designed and constructed this unit so that when used properly and for the intended purpose it does not pose any danger to people or property. Nevertheless, residual risks can remain. You should therefore observe the following notes.

### 2.1 Intended purpose

This device is a piezo-operated ultrasonic device for use in dental applications. It is mainly used for the treatment of periodontal defects. In addition, the device is used in the area of prophylaxis, peri-implantitis treatment as well as dental hygiene.

### 2.2 Intended use

The ultrasonic device is designed for use in periodontology, for the removal of plaque and for the cleaning of tooth surfaces. This is done via cavitation, polishing, grinding and scraping. Further assistance for the treatment can be provided by using hydroxyl and/or fluorapatite as a polishing agent in periodontology. Only agents approved by the manufacturer may be used. The treatment is gentle on the teeth and almost painless.

#### Paro handpiece applications

- Periodontal treatment
  - Thorough removal of biofilm and concrement, and smoothing of the root surface
- Recall
  - Removal of biofilm and gentle treatment of root surfaces even with frequent instrument usage
- Peri-implantitis treatment
  - Cleaning of implant surfaces using fibrous composite materials and special plastic instruments. Avoids damaging the implant surfaces

#### Scaler handpiece applications

- Subgingival and supragingival removal of dental calculus and concrement

The piezo-ceramic drive of the Vector Scaler allows efficient removal of deposits while providing the best and gentlest possible protection of sensitive tissue structures. The ergonomic design of the handpiece offers six long-lasting high performance LEDs to provide the best possible illumination even into areas that are difficult to see.

## 2.3 Improper use



### WARNING

#### Risk of explosion due to ignition of combustible materials

- › Do not operate the unit in any rooms in which inflammable mixtures may be present, e.g. in operating theatres.

Any other usage or usage beyond this scope is deemed to be improper. The manufacturer accepts no liability for damages resulting from this. In these cases the user/operator will bear the sole risk.

## 2.4 General safety information



### WARNING

#### Contraindication

- Ultrasound vibrations can interfere with the function of heart pacemakers and defibrillators.
- › Patients with heart pacemakers or defibrillators should not be treated with these instruments.

- › Always comply with the specifications of all guidelines, laws, and other rules and regulations applicable at the site of operation for the operation of this unit.
- › Check the function and condition of the unit prior to every use.
- › Do not convert or modify the unit.
- › Comply with the specifications of the Installation and Operating Instructions.
- › The Installation and Operating Instructions must be accessible to all operators of the unit at all times.

## 2.5 Specialist personnel

### Operation

Unit operating personnel must ensure safe and correct handling based on their training and knowledge.

- › Instruct or have every user instructed in handling the unit.

**Installation and repairs**

- › Installation, readjustments, alterations, upgrades and repairs must be carried out by Dürr Dental or by qualified personnel specifically approved and authorized by Dürr Dental.

**2.6 Electrical safety**

- › Comply with all the relevant electrical safety regulations when working on the unit.
- › Never touch the patient and unshielded plug connections on the device at the same time.
- › Replace any damaged cables or plugs immediately.

**Observe the EMC rules concerning medical devices**

- › The unit is intended for use in professional healthcare facilities (in accordance with IEC 60601-1-2). If the appliance is operated in another environment, potential effects on electromagnetic compatibility must be taken into account.
- › Do not operate the unit in the vicinity of HF surgical instruments or MRT equipment.
- › Maintain a minimum distance of at least 30 cm between the unit and other electronic devices.
- › Keep a minimum distance of 30 cm between the unit and mobile radio devices.
- › Note that cable lengths and cable extensions have effects on electromagnetic compatibility.

The following accessories can have an effect on the electromagnetic compatibility:  
 Mains cable . . . . . 9000100846  
 Flexible foot switch cable . . . . . 9000-119-130E


**NOTICE**

**Negative effects on the EMC due to non-authorized accessories**

- › Use only Dürr Dental parts or accessories specifically approved by Dürr Dental.
- › Using any other accessories may result in increased electromagnetic interference emissions or the unit having reduced electromagnetic immunity, leading to an erroneous operation mode.


**2.7 Only use original parts**

- › Only use Dürr Dental parts or accessories and special accessories specifically approved by Dürr Dental.
- › Only use only original wear parts and replacement parts.

 Dürr Dental accepts no liability for damages or injury resulting from the use of non-approved accessories or optional accessories, or from the use of non-original wear parts or replacement parts. The use of non-approved accessories, optional accessories or non-genuine wear parts / replacement parts (e.g. mains cable) can have a negative effect in terms of electrical safety and EMC.

**2.8 Transport**

The original packaging provides optimum protection for the unit during transport. If required, original packaging for the unit can be ordered from Dürr Dental.

 Dürr Dental will not accept any responsibility or liability for damage occurring during transport due to the use of incorrect packaging, even where the unit is still under guarantee.

- › Only transport the unit in its original packaging.
- › Keep the packing materials out of the reach of children.



## 2.9 Disposal

### Unit



The unit must be disposed of properly. Within the European Union, the unit must be disposed of in accordance with EU Directive 2012/19/EU (WEEE).

- › If you have any questions about the correct disposal of parts, please contact your dental trade supplier.






The unit may be contaminated. Instruct the company disposing of the waste to take the relevant safety precautions.

- › Decontaminate potentially contaminated parts before disposing of them.
- › Uncontaminated parts (e.g. electronics, plastic and metal parts etc.) should be disposed of in accordance with the local waste disposal regulations.
- › If you have any questions about the correct disposal of parts, please contact your dental trade supplier.

### 3 Overview

#### 3.1 Vector Paro / Vector Paro Pro

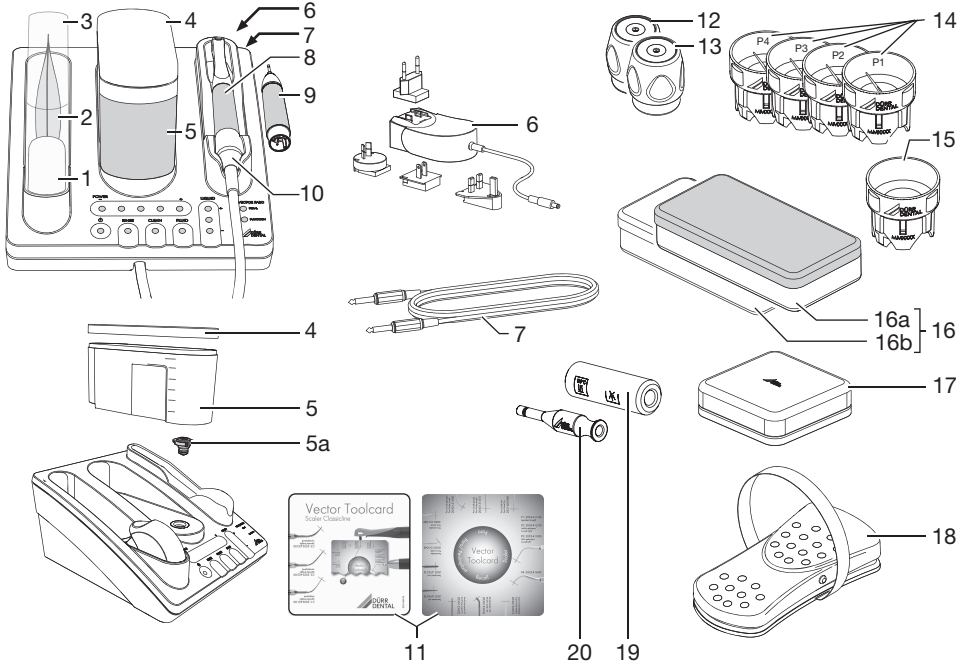

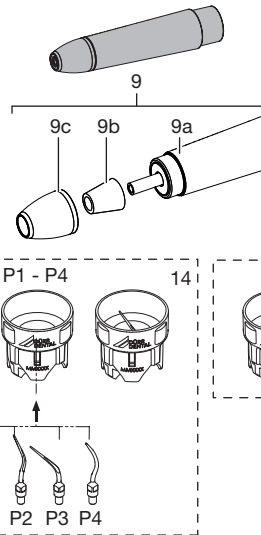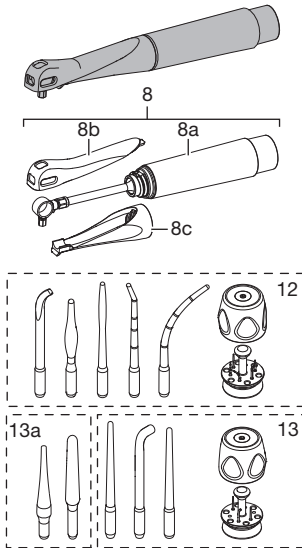

- 1 Piercing mechanism for the fluid bag
- 2 Fluid-bag "Vector Fluid Polish"
- 3 Cover for fluid bag
- 4 Cover for fluid container
- 5 Fluid container
- 5a Valve
- 6 Power supply unit
- 7 Flexible foot switch cable
- 8 Paro handpiece
- 8a Handpiece
- 8b Ring cover
- 8c Rotary adaptor
- 9 Scaler handpiece
- 9a Handpiece
- 9b Light conductor
- 9c Cover
- 10 Handpiece hose
- 11 Vector toolcard

- 12 Paro tool kit (blue ring)
- 13 Recall/Implant tool kit (black ring)
- 13a Instruments for peri-implantitis treatment, peri-implant soft and peri-implant hard
- 14 Scaler tool kit torque wrench with integrated PREMIUMLINE instrument P1 - P4
- 15 Torque wrench for PREMIUMLINE instruments (P1 - P4)
- 16 Steri-box
- 16a Paro steri-box
- 16b Scaler steri-box
- 17 Service kit
- 18 Flexible foot switch
- 19 Rinsing adapter for scaler instruments
- 20 Rinsing adapter for handpieces (yellow)

### 3.2 Scope of delivery

The following items are included in the scope of delivery (possible variations due to country-specific requirements and/or import regulations):

**Vector Paro . . . . . 2031-50**

- Paro handpiece
- Power supply unit
- Flexible foot switch (including cable)
- 3 V lithium CR 2032 battery
- Paro tool kit with integrated instruments
- Recall/Implant tool kit with integrated instruments
- Paro steri-box
- Combination wrench
- Vector toolcard
- Service kit
- Vector fluid polish MORE EFFECTIVE
- Vector cleaner, special cleaner
- Vector/RinsEndo disinfection, first application, 120 ml
- Vector Paro/Vector Paro Pro installation and operating instructions
- Quick start instructions
- Vector DVD: "Clinical Application" and "Tips and Tricks"

or

Vector Paro Pro . . . . . 2031-51

like Vector Paro, 2031-50 plus additionally

- Scaler handpiece
- Scaler instrument P1
- Scaler steri-box

### 3.3 Accessories

The following items are required for operation of the device, depending on the application:

- Paro handpiece . . . . . 2031-700-00
- Scaler handpiece . . . . . 2032-200-00
- Flexible foot switch . . . . . 2031-600-00
- Paro steri-box (cover: silver) . . . . . 2031-330-00
- Scaler steri-box (cover: blue) . . . . . 2032-330-00
- Service kit . . . . . 2031-340-00
- Vector toolcard . . . . . 2031-400-01
- Scaler instrument torque wrench . . . . . 2032100004
- Cover, Paro or Recall/Implant tool kit . . . . . 2030-150-06E
- Combination wrench . . . . . 2030-137-01E

- Rinsing adapter set (rinsing adapter for scaler instruments and rinsing adapter for handpieces (yellow)) . . . . . 2032100008

### Instruments for Paro handpiece

#### Paro tool kit

- Paro tool kit, complete . . . . . 2031-450-00
- Paro curette (3 pcs.) . . . . . 2030-151-04E
- Paro lancet (3 pcs.) . . . . . 2030-151-02E
- Paro probe plus (3 pcs.) . . . . . 2031-400-06E
- Paro probe straight (3 pcs.) . . . . . 2030-151-01E
- Paro probe curved (3 pcs.) . . . . . 2030-151-03E

#### Recall/Implant tool kit

- Recall/Implant tool kit, complete . . . . . 2031-460-00
- Recall probe straight CFRP . . . . . 2030-153-02E
- Recall curette CFRP (3 pcs.) . . . . . 2030-153-05E
- Supra probe flexible (3 pcs.) . . . . . 2030-152-01E

### Instruments for peri-implantitis treatment

- Periimplant soft (3 pcs.) . . . . . 2031-474-01E
- Periimplant hard (3 pcs.) . . . . . 2031-473-01E

### Instruments for Scaler handpiece

#### PREMIUMLINE

- Tool-Kit Scaler P1, straight . . . . . 2032-411-00
- Tool-Kit Scaler P2, right curved . . . . . 2032-412-00
- Tool-Kit Scaler P3, left curved . . . . . 2032-413-00
- Tool-Kit Scaler P4, supra . . . . . 2032-414-00

### 3.4 Consumables

The following materials are consumed during operation of the device and must be ordered separately:

- Vector Fluid Polish MORE EFFECTIVE, 200 ml . . . . . CWZ510C2350
- Vector/RinsEndo Disinfection . . . . . CDZ501C6150
- Vector cleaner, special cleaner for hose system, 4 x 2.5 l . . . . . CCA531A6150
- ID 213
- Instrument disinfection . . . . . CDI213C6150
- FD 322
- rapid surface disinfection . . . . . CDF322C6150
- FD 350 Classic
- disinfection wipes . . . . . CDF35CA0140
- FD 370 cleaner - for medical practices . . . . . CCF370C6150
- FD 366 quick-acting disinfectant for sensitive surfaces . . . . . CDF366C6150

### 3.5 Wear parts and replacement parts

The following working parts need to be changed at regular intervals (refer to the "Maintenance" section):

Instruments for Paro handpiece and Scaler handpiece, see "3.3 Accessories"

Service kit . . . . . 2031-340-00

Light conductor for Scaler handpiece (4 pcs.) . . . . . 2032-200-03E



Information about replacement parts is available from the portal for authorised specialist dealers at:  
[www.duerrdental.net](http://www.duerrdental.net).

## 4 Technical data

### Electrical data – power supply unit

Rated voltage	V AC	100 - 240
Mains frequency	Hz	50 - 60
Rated power	W	30
Current consumption	A	1 - 0.5
Protection class		II
Type of protection		IP 20

### Electrical data – basic unit and handpieces

Voltage	V DC	24
Electrical power <i>ON</i>	W	4.8
Electrical power <i>Active</i>	W	9.6
Electrical power <i>Standby</i>	W	1.2
Type of protection		IP 20

### Classification

Medical Devices Directive (93/42/EEC)	Class IIa
---------------------------------------	-----------

### General technical data – basic unit and handpieces

Operating frequency, Paro handpiece	kHz	approx. 23
Operating frequency, Scaler handpiece	kHz	approx. 27–32
Duty cycle	%	100
Fluid Polish bag contents	ml	200
Fluid container fill capacity	ml	600
Paro handpiece water consumption	ml/min	approx. 3.3
Scaler handpiece water consumption	ml/min	approx. 30–45
Max. surface temperature of instruments	°C	58

### Weight

Paro basic unit	kg	1.5
Scaler basic unit	kg	1.43
Paro handpiece	g	approx. 59
Scaler handpiece	g	approx. 56


### Dimensions (W x H x D)

Paro basic unit	cm	21.5 x 25.2 x 16
Scaler basic unit	cm	15.3 x 25.2 x 16
Paro handpiece	cm	Ø 1.9 x 14.8
Scaler handpiece	cm	Ø 2.0 x 9.4


<b>Battery for flexible foot switch</b>		
Voltage	V	3
Type		Lithium CR2032
<b>Ambient conditions during transport and storage</b>		
Temperature	°C	-15 to +60
Relative humidity	%	Max. 95
<b>Ambient conditions during operation</b>		
Temperature	°C	+10 to +40
Relative humidity	%	Max. 80
<b>Electromagnetic compatibility (EMC) Interference emission measurements</b>		
High-frequency emissions in accordance with CISPR 11		Group 1 Class B
Harmonics in acc. with IEC 61000-3-2		Not applicable
Voltage fluctuations/flickers in acc. with IEC 61000-3-3		Not applicable
<b>Electromagnetic compatibility (EMC) Interference immunity tests</b>		
Static electricity discharge in accordance with IEC 61000-4-2		Compliant
Magnetic field for a supply frequency (50/60 Hz) in accordance with IEC 61000-4-8		Compliant
Emitted HF disturbance variables in accordance with IEC 61000-4-3		Compliant

## EN 4.1 Type plate

The type plate is located on the rear of the device.

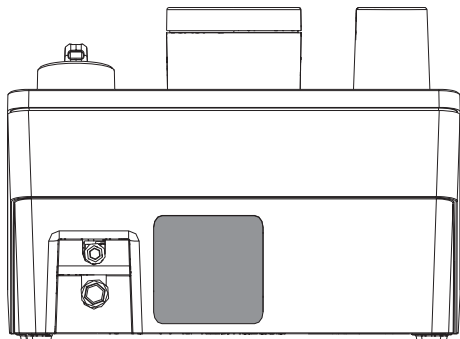

The seal is located on the base of the device.

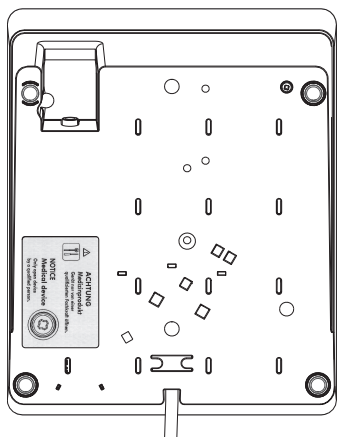



**i** Incorrectly performed work can endanger the correct operation and safety of the device.

The seal must only be removed by a qualified expert.

The device must only be opened by a qualified expert.

## 4.2 Order number and serial number for the handpieces

The serial number **SN** of the handpieces is located in the area marked in grey.

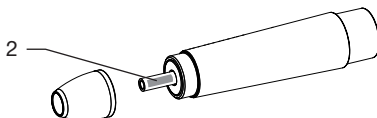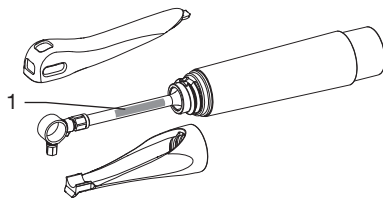

- 1 Paro handpiece
- 2 Scaler handpiece

## 4.3 ID number for tool kits

There is an ID number on the Paro and Recall/Implant tool kit covers and on the Scaler tool kits. The ID number serves to document the reprocessing.

These parts may no longer be used after a certain number of reprocessing cycles, or after the end of the service life.

The ID number is made up of the following marking: MMXXXX

- MM Date of manufacture: year and month
- XXXX Consecutive alphanumeric ID number

Paro and Recall/Implant tool kit cover



Scaler tool kit





## 4.4 Evaluation of conformity

This device has been subjected to conformity acceptance testing in accordance with the current relevant European Union guidelines. This equipment conforms to all relevant requirements.

## 5 Operation

### 5.1 Handpieces

During treatment using the Paro or Scaler handpieces, they can be operated using just water or using water and Vector Fluid Polish as required. Vector Fluid Polish can be switched on or off using the operating panel.

#### Paro handpiece

On the Paro handpiece, the Vector basic unit produces a largely linear vibration of the instrument (perpendicular to the longitudinal axis of the handpiece) with an adjustable amplitude of approx. 15–40  $\mu\text{m}$  and a frequency of approx. 23 kHz.

During treatment using the Paro handpiece fluid is applied as a pulsating jet. Once the flexible foot switch is released a small amount of fluid will still drip out. The amount of water that comes out is pre-defined automatically and cannot be changed.

#### Scaler handpiece

On the Scaler LED handpiece, the Vector basic unit generates a spatial vibration of the instrument tip (transverse to the instrument axis) with an amplitude of approx. 20–120  $\mu\text{m}$ .

During treatment using the Scaler handpiece, fluid is applied as a constant stream. The amount of water that comes out can be adjusted via the operating panel.

6 LEDs are integrated into the front section of the handpiece. As soon as the flexible foot switch is switched on and "Power" is selected on the operating panel, the LEDs light up. The LEDs go out approx. 4 seconds after the flexible foot switch is released.

#### Handpiece detection

The last used settings in the operating panel are saved by the relevant handpiece.

After the device is switched off and back on again the handpiece detection is maintained. Unplugging the mains plug of the device will cause the last used settings to be lost.

### 5.2 Fluid container

The fluid level is monitored via a sensor. When the level drops to a minimum level, the LEDs in the area of the fluid container flash and a warning signal sounds (3x).

### **5.3 Instrument change**

A torque wrench is required for inserting/changing the instruments. This is integrated in the tool kit cover for the instruments of the Paro hand-piece. For the instruments of the Scaler hand-piece there is a separate torque wrench.

## 6 Requirements

The room chosen for installation must satisfy the following requirements:

- Closed, dry room
- Clean, level and sufficiently stable subsurface
- There should be no large interference fields present (e.g. strong magnetic fields) that could interfere with correct operation of the unit.
- The required ambient conditions are satisfied (refer to the "Technical Data").


## 7 Installation

### 7.1 Establishing the electrical connections


The connection ports are located in the recess on the rear of the unit.

**Requirements:**

- ✓ Correctly installed power outlet in the vicinity of the unit (max. length of mains cable 3 m)
- ✓ The plug connection of the power supply is freely accessible so that it can be quickly disconnected in the event of danger.
- ✓ Mains voltage must match the information shown on the type plate of the power supply unit
- ✓ The supply voltage of the power supply unit matches the data on the type plate.
- › Plug in the connecting plug of the connecting cable into the connecting socket of the unit.

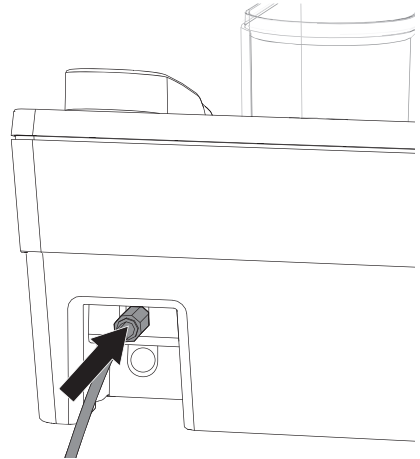

- › Fit the correct country-specific adapter.

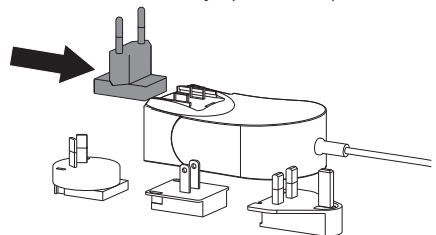

- › Plug the mains plug into the power outlet.

## 7.2 Connecting the flexible foot switch

**i** The flexible foot switch can be operated using a foot switch cable or wirelessly.


LED PEDAL flashes in orange after the unit is switched on:

- No cable connection between the unit and the foot switch.
- Pairing not performed for wireless operation.


The LED will keep flashing until a cable connection is established or pairing is performed.

### Operation with a cable

- › Plug the connector of the foot switch cable into the connecting socket on the flexible foot switch.

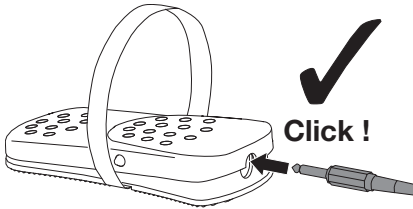

- › Plug the connector of the foot switch cable into the connecting socket on the unit.

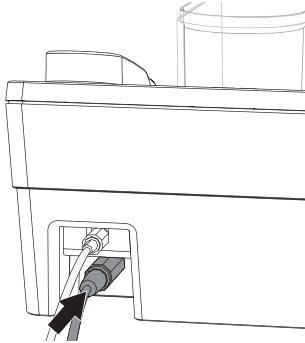

### Wireless operation

If the flexible foot switch is to be operated wirelessly, pairing (synchronization/coupling) of the flexible foot switch and the unit needs to be performed during initial start-up.

**i** In order to avoid interference in wireless operation, we recommend that a maximum of 4 flexible foot switches is used wirelessly within a single surgery. If interference does occur during wireless operation, we recommend using the flexible foot switch with foot switch cables. Wireless operation is not possible if the foot switch cable is connected to either the unit or the flexible foot switch.

### **!** CAUTION Risk of injury

Mixing up flexible foot switches can cause malfunctions such as the accidental operation of a different hand-piece. This can cause injuries.


- › If multiple devices are used simultaneously in wireless operation, make sure that the flexible foot switch that is paired with the relevant device is always used.
- › Make sure they are stored together as well.



### Performing pairing

- ✓ Place the flexible foot switch ready.
- ✓ Insert a battery in the flexible foot switch, "15.7 Inserting or changing the battery in the flexible foot switch".
- ✓ Disconnect the foot pedal cable from the unit and foot switch (if previously connected).

**i** Before carrying out the pairing process, make sure that no other Vector device with a flexible foot switch is running during the process within a range of around 10 m. Otherwise the connection may be established incorrectly.

Performing pairing between the flexible foot switch and the device:

- › Switch on the device .

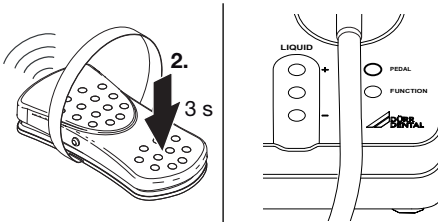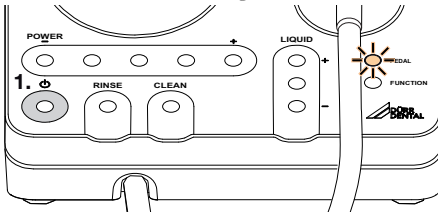

- › LED PEDAL flashes in orange.
- › Press the flexible foot switch for around 3 seconds until the orange PEDAL LED goes out.

**Result:**


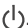
After successful pairing the unit is immediately ready for operation.



If nothing happens when the flexible foot switch is operated, it is possible that the foot switch being used is paired with a different device. In this case cancel the pairing and repeat the process.

**Canceling a pairing**

Requirements:

- ✓ No foot switch cable must be connected to the device or to the flexible foot switch.
- ✓ The orange PEDAL LED is not on or is not flashing.
- › Switch off the device .
- › Touch the LIQUID operating panel, keep touching it and switch on the device .

**Result:**

When the orange PEDAL LED flashes, this shows that the current pairing has now been cancelled.

## 8 Commissioning

### 8.1 Function check

To finish the initial start-up process, all of the connections must be checked to make sure they are securely seated and leak-tight.

Check for correct operation:

- Operating panel
- Flexible foot switch
- Visual and acoustic signals

## 9 Components

### 9.1 Paro handpiece




#### **DANGER**

**Risk of injury / treatment may not be successful**

If a handpiece is damaged by dropping it, this can cause injuries and may result in an unsuccessful treatment.


- › Perform a visual inspection for cracks and damage.
- › Any damaged parts should be replaced immediately – if necessary, replace the entire handpiece.

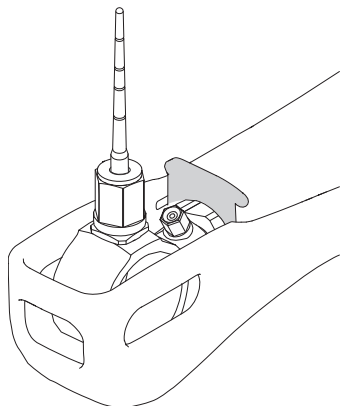

There is a sloping plastic piece on the rotary adaptor. This plastic piece is designed to hold soft tissue slightly away from the instrument, e.g. when carrying out buccal treatment on the premolars in the lower jawbone. This prevents displacement of fluid caused by the movement of soft tissue.

During operation the resonating body moves. This movement can cause friction heat if it comes into direct contact with dry mucous membrane. The ring cover prevents contact between the mucous membranes (soft tissue) and the resonating body.

#### Settings

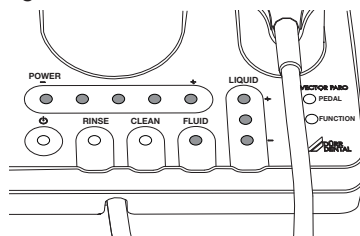

#### **POWER**

The power is adjusted in the POWER operating panel.

The Vector Paro handpiece enables optimum adjustment of the ultrasonic power according to the medical indication in conjunction with the corresponding instrument.

The power can be adjusted between 5 settings on the Scaler handpiece, "POWER".

Operating frequency:

The operating frequency of the Scaler handpiece is 23 kHz (15–35 µm).

#### **FLUID**


Vector Fluid Polish can be switched on or off using the FLUID operating panel.

#### **LIQUID**


The amount of water that comes out is predefined on the Paro handpiece and can not be changed.

#### **Disassembly**


- › Unscrew the instrument "Inserting/changing instruments".
- › Pull off the hose connection from the Paro handpiece.

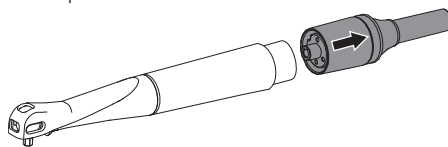

- › Press the ring cover carefully against the clip closure to detach it and then remove.

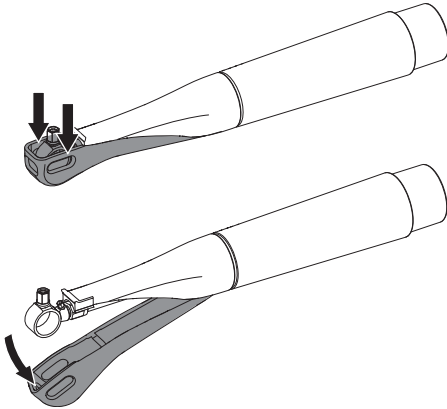


- › Twist the rotary adaptor anti-clockwise and take it off.

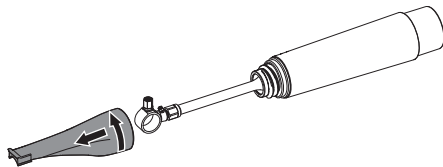

### Assembly


- › Guide the rotary adaptor over the resonating body and then twist it clockwise as far as it will go.

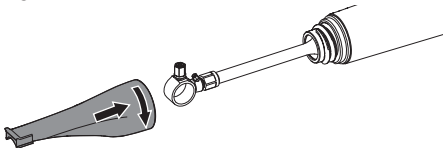


Only when the rotary adaptor has been twisted in as far as it will go is it possible to correctly insert the ring cover.

- › Insert the ring cover at an angle in the hand-piece sleeve.

- › Press the ring cover downwards until the clip closure engages.

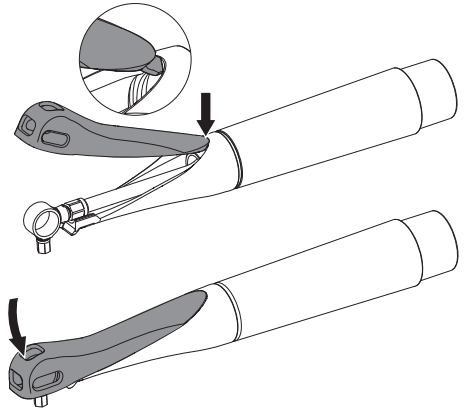

- › Connect the hose connection at the hand-piece.

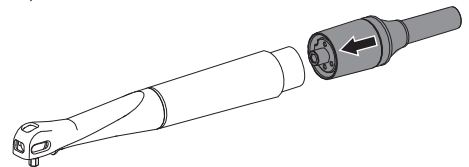

- › After assembly of the handpiece, check the nozzle for correct operation.

### Checking the operation of the nozzle





#### CAUTION

#### Risk of burns if no fluid comes out

During operation the resonating body and the instrument move. If fluid does not come out, friction heat can build up under direct contact between the resonating body or the instrument and dry mucous membranes, which can cause burns.

- › Only operate the handpiece with the nozzle inserted and working correctly.
- › Only operate the handpiece with the rotary adaptor and ring cover in proper working order.

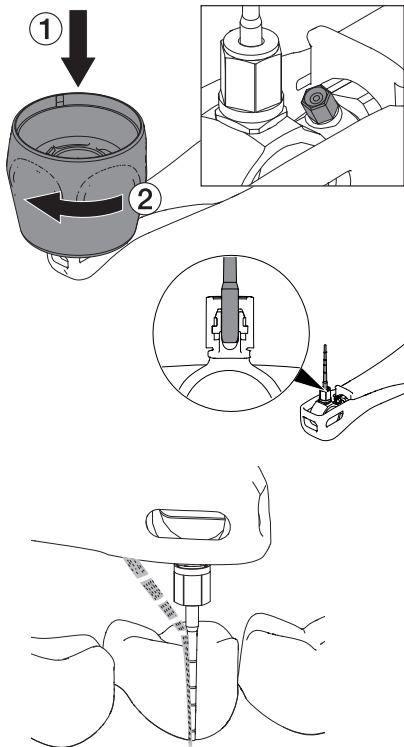

- › Preferably insert the "Paro straight probe" into the chuck and secure using the tool kit cover.
- › Switch on the device .


› Start the handpiece with the foot switch. During treatment using the Paro handpiece the nozzle ejects fluid as a pulsating stream. The fluid hits the instrument in the top third.

Possible faults:

- If spray mist is released there is possibly air in the hose lines.
- Nozzle blocked; replace the nozzle.
- Stream of fluid does not meet the instrument in the top third or completely misses the instrument; replace the nozzle.

### Instrument chuck/sealing ring check

During treatment, two sealing rings prevent fluid from entering the instrument chuck.

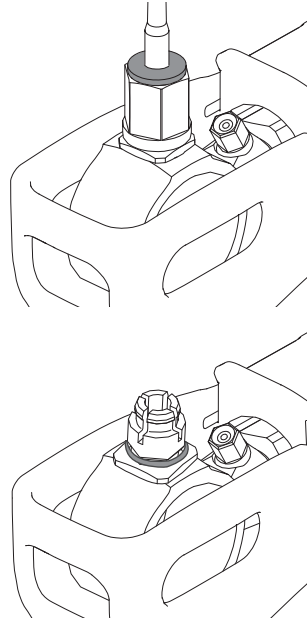

- Support seal in the union nut
- Sealing ring at the instrument chuck



- › Before every treatment check that the sealing rings are in place and that they are intact. Missing or defective sealing rings must be replaced immediately.
- › Safe and proper operation of the device can only be ensured if the instruments are correctly locked in the instrument chuck.

**CAUTION****Instrument aspiration due to missing or defective support seal**

- › Do not use the device without the support seal.
- › Check that the support seal is inserted and that it is intact.



To prevent any deformation of the instrument chuck, the union nut must only be tightened if an instrument is inserted.

- › Check that the instrument is tight and secure before placing it in the patient's mouth.

**Result:**


The following circumstances can endanger the patient and jeopardise the success of the treatment:

- Incorrectly inserted instruments
- Bent instruments
- Defective resonating ring
- Defective union nut
- Manipulation of the instrument chuck

## 9.2 Scaler handpiece

**Overview**



The Vector Scaler handpiece is suitable for the efficient removal of dental calculus and of concrement.

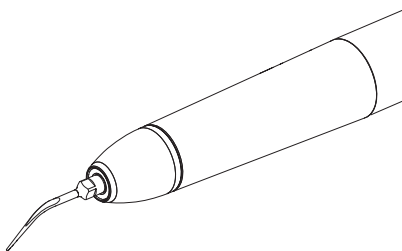

The cooling channel runs to a point just before the instrument tip. This offers the following advantages:

- Less water used, therefore less aerosol formation.
- Less contamination.
- Better clarity.
- Easier aspiration.
- Better cooling, as the fluid flows directly over the working tip.

During treatment using the Scaler handpiece, fluid is applied as a constant stream.

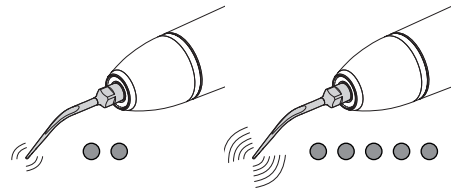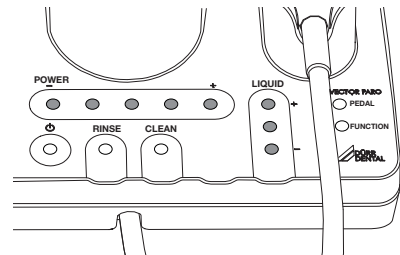**Settings****POWER**

The power is adjusted in the POWER operating panel.

The Vector Scaler handpiece enables optimum adjustment of the ultrasonic power according to the medical indication in conjunction with the corresponding instrument.

The power can be adjusted between 5 settings on the Scaler handpiece, "POWER".

Operating frequency:

The operating frequency of the Scaler handpiece is in the range between 27–32 kHz (20–120 µm)

**LIQUID**

On the Scaler handpiece the water amount can be selected from 3 settings:


LED display	Water amount
1	30 ml/min
2	37–40 ml/min

**LED display**      **Water amount**

3                      45 ml/min

**Illumination**


**Power**

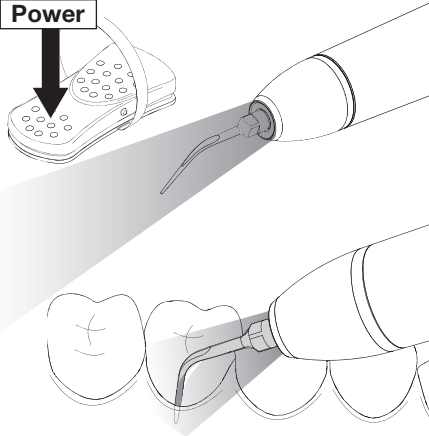

6 LEDs are integrated into the front section of the handpiece. The light conductor is located under the front cover. As soon as the flexible foot switch is operated and "Power" is selected on the operating panel, the LEDs light up. The LEDs go out approx. 4 seconds after the flexible foot switch is released.


The LEDs are actuated individually so that if one LED stops working the light source remains lit.

**Disassembly**

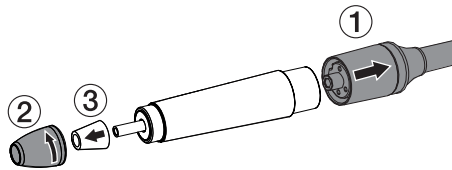

- › Dismantle the instrument.
- › Pull off the hose connection from the handpiece.
- › Unscrew the front cover from the handpiece working anti-clockwise.
- › Disconnect the light conductor.

**Assembly**

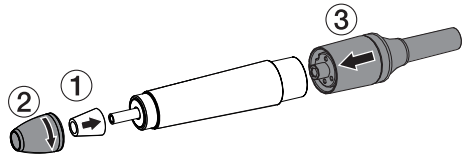


- › Connect the light conductor.
- › Screw the front cover onto the handpiece working clockwise.
- › Connect the hose connection at the handpiece.

**9.3 Instruments and tool kits**


**Overview**

Instruments of various shapes, lengths and materials are available. These are grouped according to their different applications and arranged in the tool kits.

**i** These instruments are specially designed for use with the Vector unit. No other instruments must be used.

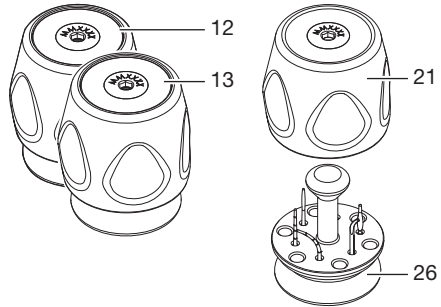

- 12 Paro tool kit
- 13 Recall/Implant tool kit
- 21 Tool kit cover
- 26 Instrument tray

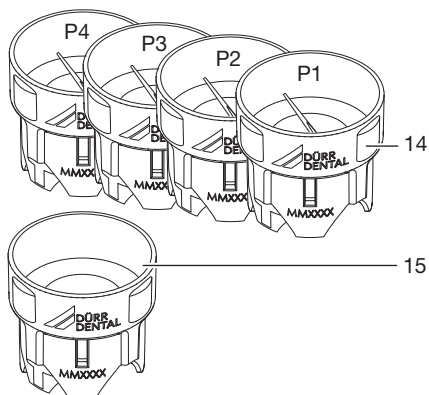

- 14 Scaler tool kit torque wrench with integrated PREMIUMLINE instrument  
 15 Torque wrench for all PREMIUMLINE instruments (P1 - P4)

The tool kits are designed for the storage, cleaning, disinfection and sterilisation of the instruments.

#### Paro and Recall/Implant tool kit

The instruments are placed in the instrument tray and sealed with the tool kit cover. The tool kit cover serves also as a torque wrench for changing the instrument.

#### Scaler tool kit

Each instrument is kept in its own separate tool kit. The tool kit cover serves as a torque wrench for changing the instrument.

#### Instruments made of metal

The use of metal instruments enables a higher application of energy.

Application areas:

- Periodontal initial treatment
- Removal of concrement and dental tartar

#### Instruments made of fibre composite (CFRP)


Instruments made of CFRP are semi-flexible.

Application areas:

- Supporting periodontitis treatment
- Gentle removal of subgingival biofilm and supragingival plaque
- Removal of discoloration under maximum protection of delicate structures, such as root cement, exposed dentine surfaces or bone islands, as well as the sensitive surfaces of teeth, dentures or implants.
- Use on non-metallic dentures

#### Paro tool kit


The Paro tool kit contains instruments for the Vector Paro handpiece for initial periodontal treatment.

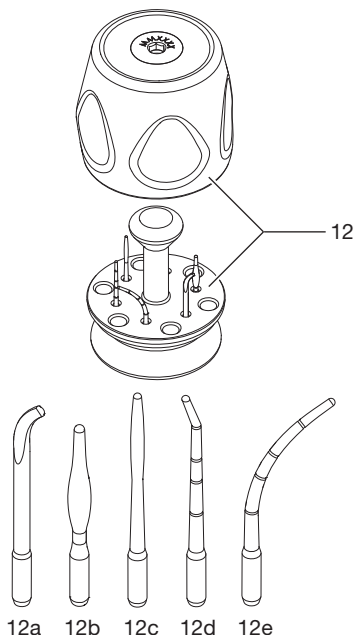

- 12 Paro tool kit  
 12a Paro curette  
 12b Paro lancet  
 12c Paro probe plus  
 12d Paro probe straight  
 12e Paro probe curved

**Recall/Implant tool kit**

The Recall/Implant tool kit contains instruments for the Vector Paro handpiece.

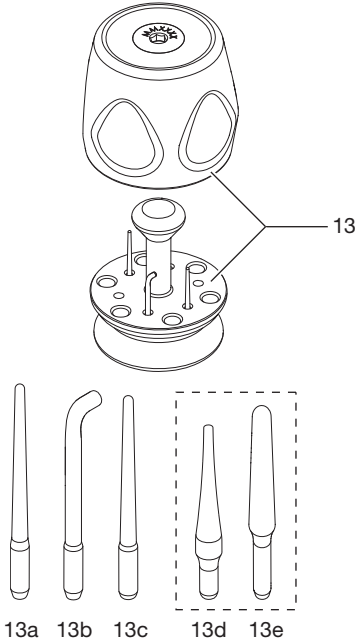

**Recall**

- 13 Recall tool kit
- 13a Recall probe straight CFRP
- 13b Recall curette CFRP
- 13c Supra probe flexible

**Implant**

Instruments for peri-implantitis treatment for sensitive implant surfaces

- 13d Periimplant soft
- 13e Periimplant hard

**Scaler tool kit**


The use of metal instruments enables a higher application of energy.

Application areas:

- Periodontal initial treatment
- Removal of concrement and dental tartar

**PREMIUMLINE instruments**

Each instrument is contained in its own tool kit. The tool kit cover serves as a torque wrench for changing the instrument.

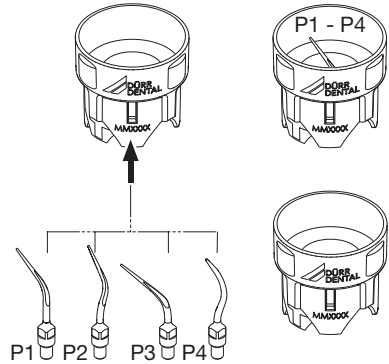

- **Scaler tool kit P1**

30 µm, straight, for the removal of subgingival deposits with a pocket depth of up to 4 mm

- **Scaler tool kit P2**

60 µm curved right, for the removal of subgingival deposits

- **Scaler tool kit P3**


60 µm curved left, for the removal of subgingival deposits

- **Scaler tool kit P4**


120 µm, for supragingival removal of coatings on smooth surfaces and for interdental areas

## 9.4 Steri-box


**Paro steri-box** (cover: silver)

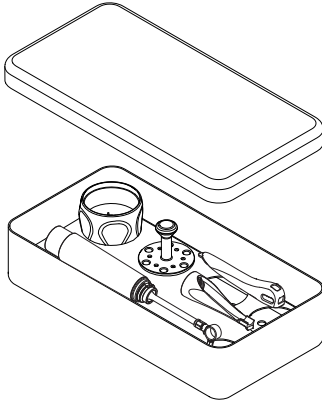

**Scaler steri-box** (cover: blue)

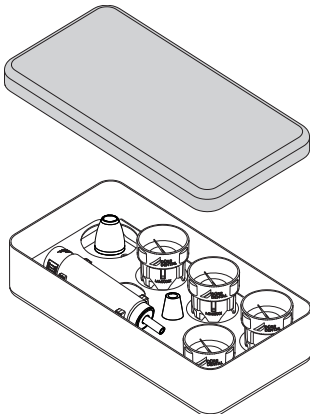


All sterilisable parts of the Vector device can be optimally placed into the steri-box.  
For steam sterilisation, the steri-boxes are placed in the small steam steriliser, "13.10 Steam sterilising".  
If sterile storage is required, wrap the steri-box in suitable sterilisation packaging in accordance with DIN11607-1 and seal.

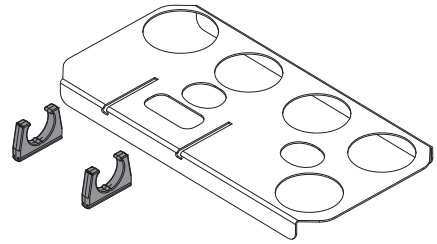

The rubber clips on the carrier plate can be replaced if necessary.

## 9.5 Flexible foot switch

The handpieces are operated with the flexible foot switch.

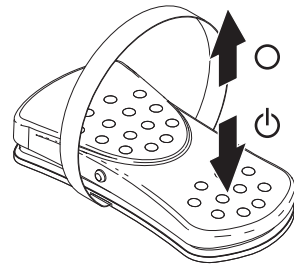

In wireless operation the flexible foot switch is supplied with voltage from a battery. If the battery power starts to fail, the orange PEDAL LED on the unit lights up.




Interference can occur in wireless operation if the battery power is low, so it is important to change the batteries in good time.

The service life of the battery is around 1 year or approx. 900 treatments.

If the battery is empty or not present then the flexible foot switch can be connected with a cable to the unit so that the treatment can be continued. The battery can then be inserted later on.


Changing the battery – "15.7 Inserting or changing the battery in the flexible foot switch".

## 9.6 Service kit

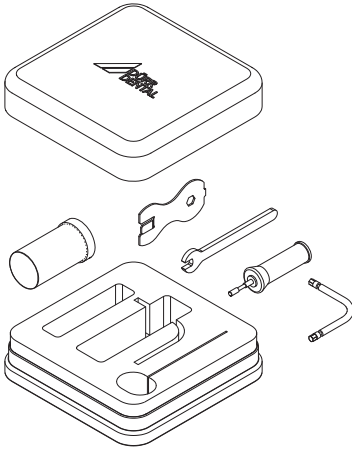

The service kit has been configured for the Vector System. The service kit contains replacement parts and tools that can be used for maintenance and repairs, "15 Maintenance".

## 9.7 Fluid container

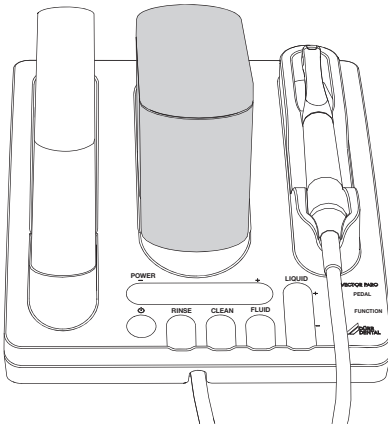

The fluid container can hold up to approx. 600ml of water or watery fluids.

- › The fluid container should be emptied and cleaned every evening and descaled as required, "12.3 Cleaning the fluid container".

**i** The fluid container is not suitable for thermal disinfection or sterilisation.

## Water quality

The water quality must meet the general requirements for water supply systems in a dental surgery and all applicable standards.

## Active ingredient solutions



### NOTICE

**Blockage of the nozzle when mixing water-based active ingredient solutions with Vector Fluid Polish.**

- This combination leads to flocculation of the fluid and to blockages as a result.
- › Vector Fluid Polish must not be added when using water-based active ingredient solutions.
- › After using a water-based active ingredient the unit should be thoroughly rinsed using warm water.

In order to widen the range of possible treatments, additional active ingredients can be added to the water.

Possible active ingredients include e.g. chlorhexidine digluconate solutions up to a total concentration of 0.2 %; other solutions may damage the system and can lead to unsuccessful treatment.

## 9.8 Vector Fluid Polish



**CAUTION**

**Extremely sensitive patients may experience reactions on their mucous membranes.**

Vector Fluid Polish contains PHB ester as a preserving agent.

- › In cases of known patient hypersensitivity to PHB esters, treatment should be carried out using just water or a water-based active ingredient solution.

Vector Fluid Polish is a polishing fluid with hydroxylapatite particles (average size <math>< 10 \mu\text{m}</math>) for the following applications:

- Smoothing work with a polishing effect
- Thorough removal of subgingival, adherent biofilms. The energy efficiency of the Vector Paro is improved by the particles added to the polish.
- Teeth cleaning
- Peri-implant mucositis and peri-implantitis
- Applications on the periodontium

Vector Fluid Polish is a ready-to-use solution, filled in a bag. Special plastic inserts on the inside of the fluid bag ensure complete emptying.



Store Vector Fluid Polish at room temperature.

To prevent it from drying out, do not expose it to direct sunlight.

Do not store it in a refrigerator. Cold Fluid Polish is highly viscous and uncomfortable for the patient.

## 9.9 Cleaning components

### Vector/RinsEndo Disinfection


Ready-to-use aldehyde-free solutions for fast-acting disinfection and cleaning of the Vector hose system. The undiluted solution is poured into the fluid container. These steps must be carried out every time a new fluid bag is inserted and before breaks in treatment of more than 24 hours, "14 Treatment breaks for more than 24 hours".

### Vector cleaner

Ready-to-use solution for the removal of acid-soluble deposits in the hose system and in the handpieces of the Vector system. Special cleaner with intensive cleaning effect and very good material compatibility.

## 9.10 Vector toolcard

The instruments are subject to wear which will depend on their material, the surfaces being treated and length of time in use.

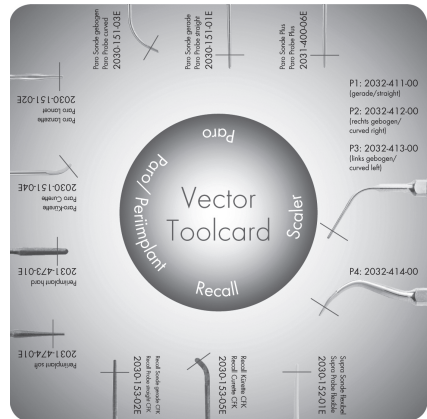

All instruments must be regularly checked with the Vector toolcard in order to check the degree of wear, "15.4 Checking instrument wear".

# EN 10 Operation

## 10.1 Display/handling



**Paro handpiece** inserted:  
 POWER adjustment possible  
 LIQUID adjustment not possible

**Scaler handpiece** inserted:  
 POWER and LIQUID adjustment possible

**No handpiece** inserted:  
 The settings cannot be changed.  
 The LEDs are only active and light up when a handpiece is inserted.  
 The last settings set up for the handpiece that is inserted remain active after the unit is switched off and back on again.  
 Any running cleaning or rinsing processes can be stopped by touching the corresponding button again.









- LED off
- LED on
- LED flashing

<b>40</b>	<b>ON / Standby</b> For switching on or off; touch the button for at least 2 seconds.
	Standby, unit switched off.
	ON, unit switched on.

	If no function is used for a period of 30 minutes, the unit will automatically switch off (standby).
<b>41</b>	<b>POWER</b> (power setting)
	1 – 5 LEDs light up depending on the power setting selected (5 LEDs = maximum power) <b>Tip:</b> The power setting can also be adjusted during treatment.
<b>42</b>	<b>RINSE</b> (rinsing/disinfecting)
	To start the rinsing process: touch the button for at least 2 seconds. The blue LED flashes during the rinsing process.
<b>43</b>	<b>CLEAN</b>
	After approx. 30 operating hours the blue LED will light up continuously – this shows that it is time for cleaning.
	To start the cleaning process: touch the button for at least 2 seconds. The blue LED flashes during the cleaning process. <b>Recommendation:</b> Clean the unit every four weeks or, at the latest, as soon as the LED lights up continuously.
<b>44 / 49</b>	<b>FLUID</b>
	To switch the fluid on or off: touch the button for at least 1 second.
	During fluid delivery the LED in the piercing mechanism also lights up.
	To activate fluid pre-delivery to the handpiece: touch the button for at least 2 seconds.
	During fluid pre-delivery the LED in the piercing mechanism also lights up.
<b>45</b>	<b>LIQUID</b> This display is only active when the Scaler handpiece is inserted.
	1 LED on = minimum fluid consumption (30 ml/minute) 3 LEDs on = maximum fluid consumption (45 ml/minute)
<b>46</b>	<b>PEDAL</b> (flexible foot switch)

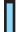


	LED on: low battery power – change the battery of the flexible foot switch.
	LED flashing: no flexible foot switch connected (cable operation) or paired (wireless operation).
47	<b>FUNCTION</b>
	LED flashing: treatment was stopped. Clean the instrument chuck with the air and water syringe and dry it, then continue the treatment.
	LED flashing: the vibration behaviour of the instrument is impaired. <ul style="list-style-type: none"> <li>– Contact pressure of the instrument is too high during treatment – reduce the contact pressure.</li> <li>– Check the instrument for signs of wear or bending.</li> <li>– Clean the Paro handpiece instrument chuck and resonant ring using the air and water syringe and dry.</li> </ul>
48	<b>LED in the fluid container</b>
	LED on: normal operation, fluid container filled enough.
	LED flashing: fluid level low. If the minimum fluid level is reached, the LED of the fluid container will start to flash and an acoustic warning signal will sound (3x high-pitched sound).

**Acoustic signals** 

Audible signals	Trigger/situation
Clicking noise	<ul style="list-style-type: none"> <li>– Touching of the operating panel</li> <li>– Function activated, e.g. RINSE, CLEAN, FLUID</li> </ul>
Long, low audible signal	– Function cannot be performed
Warning signal, 3x high-pitched audible signal	<ul style="list-style-type: none"> <li>– Fluid level at minimum</li> <li>– Waiting for fluid during the cleaning process</li> </ul>

## 10.2 Adjustment options

### Operating panel

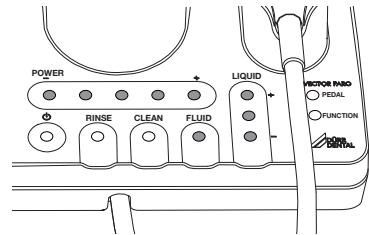

### Settings on operating panel


The settings are changed by touching but without applying pressure.



The LEDs for POWER and LIQUID are only active and light up when a handpiece is inserted.

If no handpiece is inserted, the settings for POWER and LIQUID cannot be adjusted.

### ON/Standby

The unit is switched on and switched to standby mode via the operating panel  ON / Standby. The unit is equipped with an automatic standby mode in order to save electricity. When the unit is left unused for 30 minutes it switches off automatically.

### POWER

The power is adjusted from 1 to 5; the settings are displayed via the 5 LEDs:

LED display	Power in %
1	20
2	40
3	60
4	80
5	100

The power is set to level 5 as the default factory setting on delivery. The setting for the operating frequency depends on the inserted handpiece.  
 Paro handpiece: 15–35  $\mu\text{m}$   
 Scaler handpiece: 20–120  $\mu\text{m}$

**LIQUID**

On the Scaler handpiece the water amount can be selected from 3 settings, which are displayed via 3 LEDs:

LED display	Water amount in ml/min
1	approx. 30
2	approx. 37
3	approx. 45

On the Paro handpiece the water amount is fixed and cannot be adjusted.

**RINSE**

After every treatment the system must be rinsed with water.

The rinsing process is started by touching the RINSE button; it ends automatically after approx. 30 seconds.

If a rinsing process is running it can be interrupted at any time by touching the RINSE button.

**CLEAN**

The cleaning process lasts around 10 minutes. During this step *Vector cleaner* cleaning fluid is pumped through the lines, whereby these are cleaned of any deposits.

Cleaning can be started whenever required.

Once the process is started it runs automatically until the program is finished.

We recommend that cleaning is performed every 4 weeks.

After an operating time of approx. 30 hours the blue LED on the operating panel lights up to show that cleaning needs to be performed.

The blue LED goes out when the cleaning process is completely finished.

If the cleaning process is not properly completed, the blue LED will light up every time the unit is switched on.

**FLUID**

A sensor in the piercing mechanism registers whether a fluid bag is inserted in the unit. If no fluid bag is inserted the unit will automatically operate with water (fluid).

**PEDAL (flexible foot switch)**

Operating the flexible foot switch will activate the handpiece.

If the orange LED lights up, the battery power must be checked, "15.7 Inserting or changing the battery in the flexible foot switch".

If the orange LED flashes then no flexible foot switch is connected or paired.

**FUNCTION**

If the orange LED lights up, the contact pressure of the instrument is too high or the handpiece needs to be checked.

**10.3 Preparing the device for treatment****Switch on the unit.****WARNING****Risk of cross contamination**

- › All parts must be reprocessed before every treatment.
- › If there has been no treatment for more than 24 hours then the complete fluid system must be disinfected.

› Switch on the unit.

**Result:**


Blue LED is on – unit is ready for operation.

**Inserting the fluid bag****CAUTION****Risk of injury due to sharp needle**

- There is a risk of injury when touching the needle of the piercing mechanism.
- › Proceed with caution when changing the fluid bag.


› Remove the cover of the fluid bag.


- › Check that the rubber seal is correctly seated in the device. The rubber seal must remain in position in the device during operation.

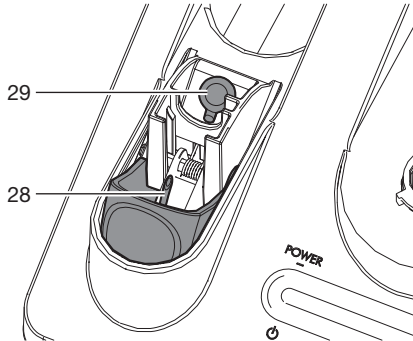


- 28 Piercing mechanism
- 29 Rubber seal

- › Before placing a new fluid bag in the unit, shake well.
- › Push the piercing mechanism to the rear.

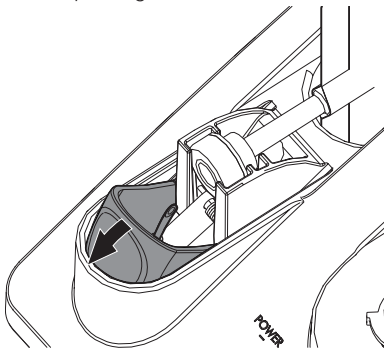


-  The closure of the fluid bag must not be removed.

- › Insert the fluid bag in the correct position in the holder so that the ball on the closure points upwards.

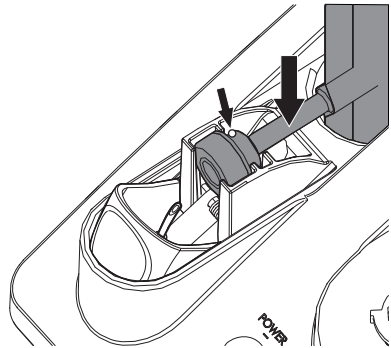

- › Slide the piercing mechanism forwards as far as it will go. In the process the needle in the piercing mechanism penetrates the closure of the fluid bag, allowing fluid to be transported via the hose to the handpiece.

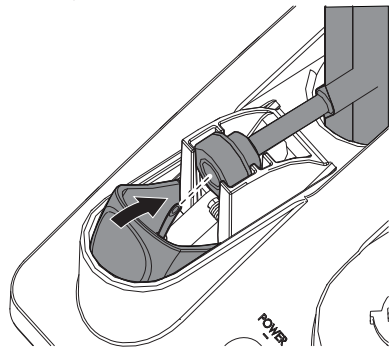

- › Fit the cover of the fluid bag.

#### Fluid pre-delivery

In the following situations the fluid line will be empty and will need to be recharged:

- ✓ After the reprocessing process
- ✓ After the unit has been out of use for an extended period of time
- ✓ After a new fluid bag has been inserted
- › Touch the Fluid button for at least 2 seconds.


#### Result:

Fluid pre-delivered to the handpiece. This step takes approx. 11 seconds and ends automati-


cally. If necessary the process can be stopped earlier by pressing the button again.

### Inserting the fluid container

- › Check whether the fluid container is full. If necessary top up the fluid container to the upper marking with lukewarm water (approx. 30 °C).
- › Insert the fluid container vertically in the correct position into the unit and press downwards until you feel it lock into position.

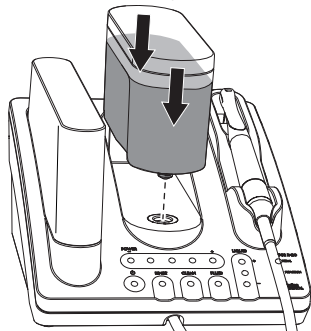

### Attaching the handpieces

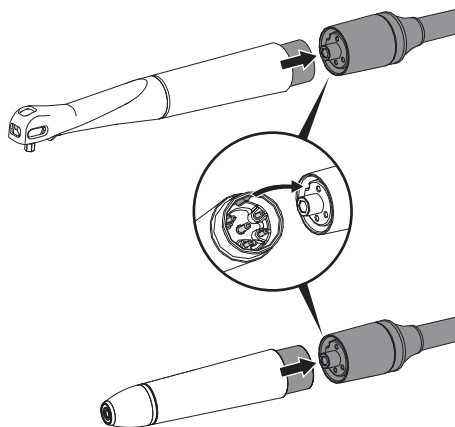

#### *Paro handpiece:*

- › Attach the Paro handpiece to the hose connection.
- › Check the ring cover for signs of damage before using it.

#### *Scaler handpiece:*

- › Attach the Scaler handpiece to the hose connection.

### Inserting/changing instruments



#### **CAUTION**

**If an attempt is made to perform a treatment with a damaged or worn instrument, this can cause injury and may result in an unsuccessful treatment.**

- › Worn or bent instruments should be replaced immediately and not reused!
- › Check the Paro handpiece instruments for signs of discolouration or roughness in the region where they are held in the chuck; if necessary dispose of them.





The tool kit cover of the instrument is used as a torque wrench for tightening the union nut of the instrument chuck.

- › Select a suitable instrument for the corresponding treatment.


#### *Instrument for Paro handpiece:*

- › Loosen the union nut approximately a 1/4 turn.
- › Insert the instrument as far as it will go into the instrument chuck.

- › Use the tool kit cover to tighten the union nut of the instrument chuck until the torque limit is reached and ratcheting of the tool kit cover can be heard.

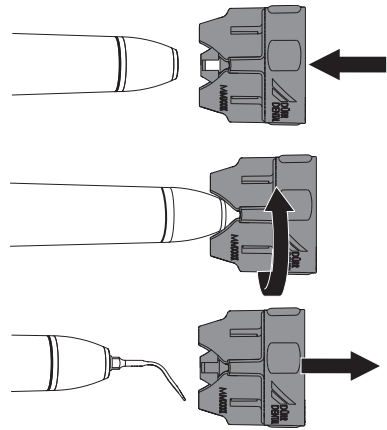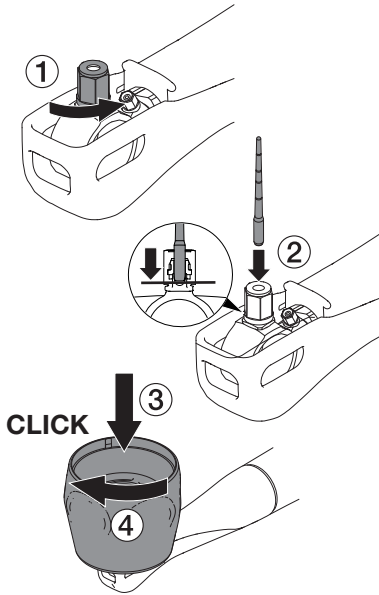

*Instrument for Scaler handpiece:*

 To prevent the instrument from being overtightened, the torque wrench slips when the required torque is reached. No ratcheting sound can be heard.

- › Always use the torque wrench to screw the instrument on and off.  
When screwing the instrument on, slowly turn the torque wrench approx. a quarter turn beyond the resistance point.

## 11 Treatment

### 11.1 Preparation

Perform the following steps before every treatment:

- › Make sure that only handpieces and instruments are used that have been reprocessed since the last treatment.
- › Check the nozzle for correct operation, "Checking the operation of the nozzle".
- › Check that the sealing ring on the instrument chuck and the support seal on the union nut are correctly seated and in perfect condition, "Instrument chuck/sealing ring check".
- › Check that the instrument is correctly seated and in perfect condition, "Inserting/changing instruments".
- › Check the fill level of the fluid container.
- › Adjust the power on the operating panel as required, "POWER".

### 11.2 Treatment using the Paro handpiece

During treatment using the Paro handpiece fluid is applied as a pulsating jet. Once the flexible foot switch is released a small amount of fluid will still drip out.




#### CAUTION

##### Risk of burns if no fluid comes out

During operation the resonating body and the instrument move. If fluid does not come out, friction heat can build up under direct contact between the resonating body or the instrument and dry mucous membranes, which can cause burns.

- › Only operate the handpiece with the nozzle inserted and working correctly.
  - › Only operate the handpiece with the rotary adaptor and ring cover in proper working order.
- › Activate the handpiece by operating the flexible foot switch.

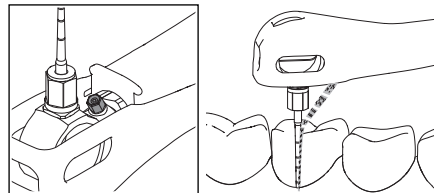




If there is too much contact pressure on the instrument during treatment, the process will be interrupted and the orange FUNCTION LED will light up.

Reduce the contact pressure and press the flexible foot switch again to continue the process and switch off the LED. If necessary, the instrument might require spraying and drying using the air and water syringe.

- › Always guide the Vector Paro instrument parallel to the root surface.


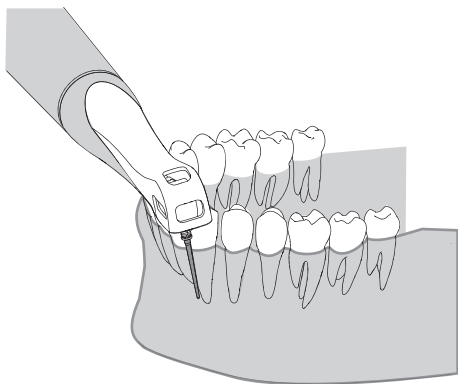


## CAUTION

### Risk of burns due to friction heat

Friction heat is created during contact between soft tissue and the union nut, which can cause burns.

- › During any treatment the soft tissue, e.g. cheeks, lips, tongue etc. must be kept well clear of the union nut.
- › Excess fluid should preferably be aspirated using the smaller saliva ejector in the dorsal oral cavity area on the contralateral side. To maintain the energy efficiency of the fluid and the Vector Fluid Polish, do not aspirate directly at the site of the treatment. If you aspirate directly at the site of treatment then no energy can be dissipated and friction heat can build up.

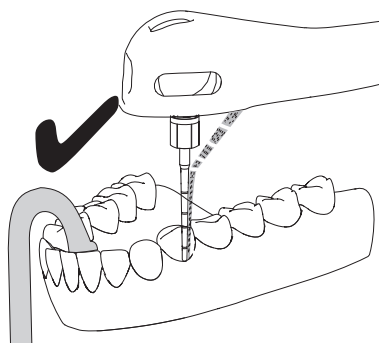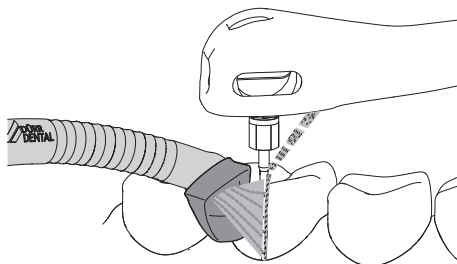

With various predecessor types of the Vector series it was possible to use Vector Fluid Abrasive.

When using Vector Paro only Vector Fluid Polish must be used, as otherwise blockages can occur.

## EN 11.3 Treatment with a Scaler hand-piece




### CAUTION


#### Risk of injury due to burns

During operation some of the components inside the handpiece heat up. Contact with hot parts can cause burns.

- › Only operate the Scaler handpiece with the cover fitted and intact.

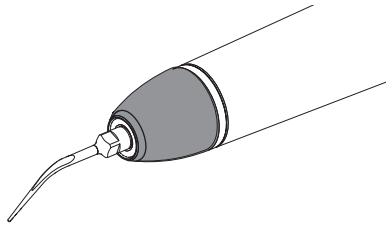

During treatment using the Scaler handpiece, it can be operated using just water (LIQUID) or using water with the addition of Vector Fluid Polish (FLUID) as required.

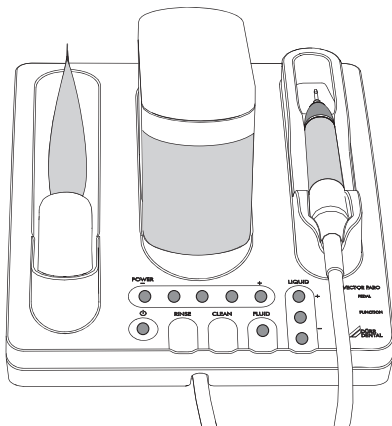

During treatment using the Scaler handpiece, fluid is applied as a constant stream.



The addition of Vector Fluid Polish reduces pain intensity for the patient.




### NOTICE

Residue of Vector Fluid Polish can cause blockages in the device.

- › Rinse the unit with water after every treatment in which Vector Fluid Polish is used (RINSE).



## 11.4 Use of the Scaler instruments


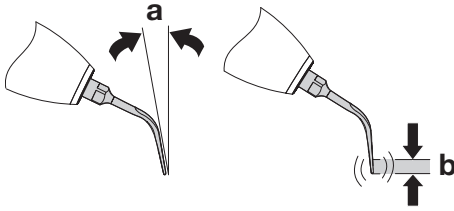


- A Contact angle approx. 10°  
 b Working area 2 mm

The active working area of the instrument is in the area of the front 2 mm.

Thanks to the minimal pain generation, treatment with the Vector Scaler can also be carried out on acute, painful parodontopathies.

- › Activate the handpiece by operating the flexible foot switch.
- › Hold the instrument at an angle of approx. 10° against the tooth and work away from the tooth.

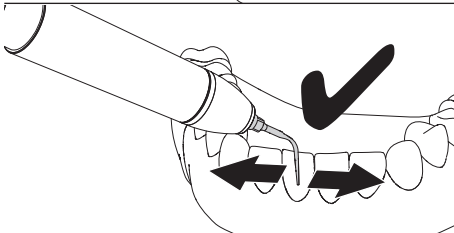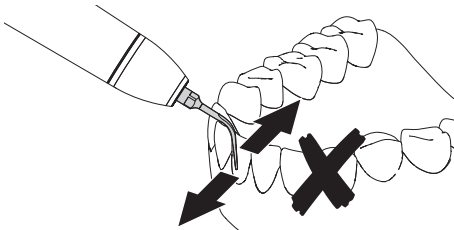

- › Keep the instrument constantly in motion: in the longitudinal direction of the tooth or transversely across the approximal surface lingually or buccally away from the tooth.
- › Guide the instrument with minimal pressure in such a way that the movements of the tip are always parallel to the tooth surface.
- › Make sure that only the side surfaces of the instruments are used. Never use the front or rear surface of the instruments.
- › Perform effective and targeted suction to provide good visibility of the field of treatment.


### Application areas

#### CAUTION Risk of injury



Accidental activation or uncontrolled activities of the handpiece can cause injuries.

- › Insert the handpiece in the handpiece holder when it is not in use.
- › Dismantle the instrument or push on the torque wrench.



- **Scaler instrument P1**, 30 µm, straight, for the removal of subgingival deposits in deep gingival pockets (up to 4 mm).  
 POWER  
 2 LEDs: 40 % operating performance

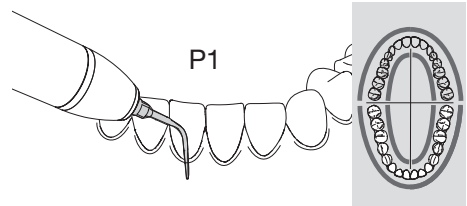


- **Scaler instrument P2**, 60 µm, curved right, for the removal of subgingival deposits  
 POWER  
 2–4 LEDs: 40 % - 80 % operating performance

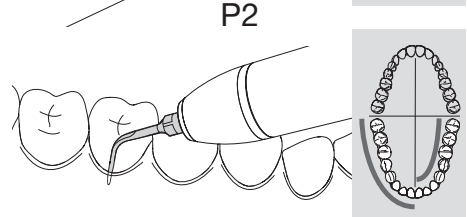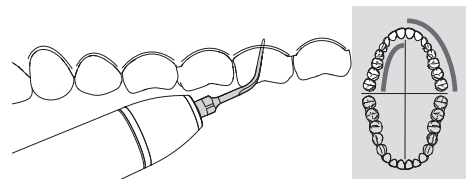

- **Scaler instrument P3**, 60 µm, curved left, for the removal of subgingival deposits  
 POWER  
 2–4 LEDs: 40 % - 80 % operating performance

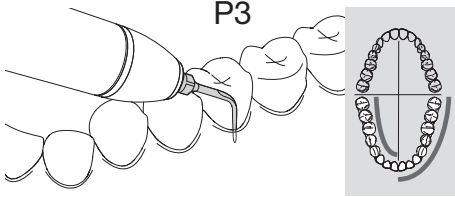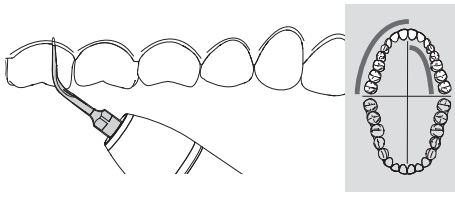

– **Scaler instrument P4**, 120 µm, for the removal of supragingival coatings on smooth surfaces and for interdental areas.  
**POWER**  
 2–3 LEDs: 40 % - 60 % operating performance

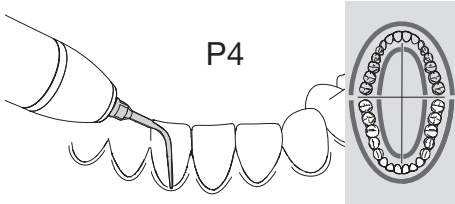

## 11.5 After every treatment

### End of treatment

- › Dismantle the instrument with the tool kit cover.
- › After every treatment, all parts used must be cleaned, disinfected and, if necessary, sterilised, "13 Reprocessing".

## 12 Cleaning

### 12.1 Cleaning of the outside surfaces

All outside surfaces must be cleaned and disinfected if they are contaminated or soiled.

- Surface of the device
- Handpiece hose
- Protective cap of the fluid polish bag
- Fluid container
- Tool kit for Scaler without instruments

As the surface disinfectant, we recommend using a disinfectant that is compatible with the materials and meets the general dental hygiene standards, e.g.:

- Dürr Dental FD 322 rapid surface disinfectant
- Dürr Dental FD 350 disinfection wipes
- Dürr Dental FD 366 sensitive rapid surface disinfectant




### NOTICE


**Liquid can cause damage to the unit.**

- › Do not spray the unit with cleaning and disinfectant agents.
- › Make sure that liquid does not get inside the unit.

- › For pre-cleaning, remove coarse organic soiling with a paper towel.
- › Clean the surfaces with a moist, soft, lint-free cloth.

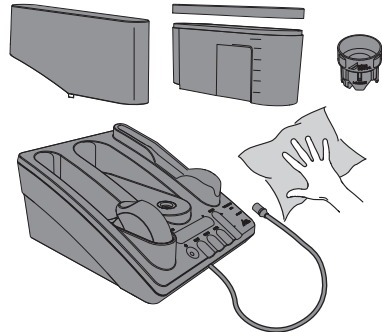


## 12.2 Activating the cleaning process of the device

 We recommend that cleaning is performed every 4 weeks. The cleaning process can be started at any time as required.


After approx. 30 operating the CLEAN LED will light up on the operating panel to show that cleaning is required.

A cleaning cycle comprises two steps that need to be started one after the other:


- ✓ CLEAN
- ✓ RINSE (rinsing/disinfecting)

 The cleaning process is not finished until both steps have been successfully performed one after the other.


### CLEAN:

 Before starting the CLEAN or RINSE (rinsing/disinfecting) functions, if there is a fluid bag in the unit it will need to be removed. If an attempt is made to start the function with a fluid bag inserted, a long, low audible signal will sound and the action will not be performed.

- › Remove the cover of the fluid bag.
- › Push the piercing mechanism to the rear.
- › Take out the fluid bag.
- › Slide the piercing mechanism forwards as far as it will go.

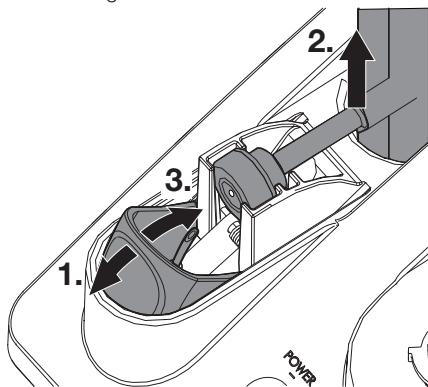


- › Pull off the hose connection from the handpiece.

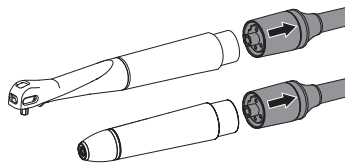


- › Place the handpiece hose in a sink or in a suitable container.



- › Pour 2 sealing caps (approx. 40 ml) of Vector cleaner undiluted into the empty fluid container.

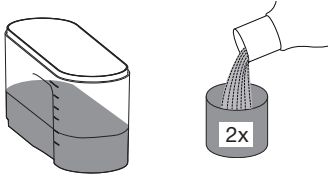


- › Touch the CLEAN button for at least 2 seconds.

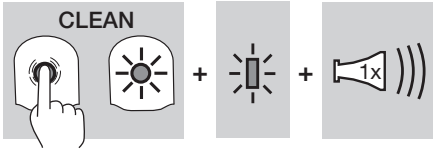

LED CLEAN and the LED in the fluid container flash and a clicking sound is heard.

The device is cleaned for approx. 10 minutes with *Vector cleaner* until the fluid container is empty. The process ends automatically.

LED RINSE lights up and a cyclically repeating warning signal sounds as an indication that the device must be rinsed with water after the cleaning with *Vector cleaner*.

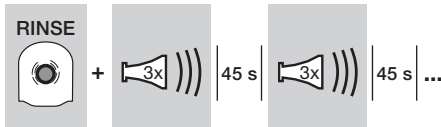

LED CLEAN and the LED in the fluid container flash and a clicking sound is heard.

The device is cleaned for approx. 10 minutes with *Vector cleaner* until the fluid container is empty. The process ends automatically.


LED RINSE lights up and a cyclically repeating warning signal sounds as an indication that the device must be rinsed with water after the cleaning with *Vector cleaner*.

**RINSE (rinsing/disinfecting):**



As a result of the rinsing with water the special cleaner *Vector cleaner* is removed from the system. Any residue of the cleaning agent could cause irritation in the patient.


- › Top up the fluid container to the upper marking with water.

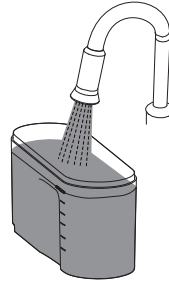

- › Touch the RINSE (rinsing/disinfecting) button for at least 2 seconds.

LED RINSE flashes.

**RINSE**

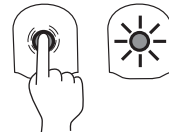

The device is rinsed with water for 30 seconds. The process ends automatically.

After the end of the complete cleaning process (CLEAN + RINSE) the CLEAN LED goes out again and an acoustic signal sounds (3x high-pitched audible signal).

If the cleaning process is not performed in full or it is stopped prematurely, then the blue CLEAN LED will come on every time the device is switched on.


## 12.3 Cleaning the fluid container

The fluid container should regularly be cleaned and descaled.

How soon descaling is required depends primarily on the hardness of the water used. However, at the latest descaling must be performed when first signs of limescale become apparent.

### Cleaning:

- › Fill up the fluid container to the upper marking with cleaning solution.

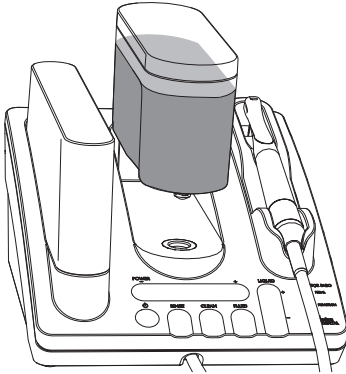


- › Allow the cleaning agent to soak as directed in the manufacturer's information.
- › Completely empty the fluid container.
- › Thoroughly rinse the fluid container with water and dry.

### Descaling:


- › Fill up the fluid container to the upper marking with e.g. 10% citric acid solution.
- › Allow the descaler to soak in; follow the manufacturer's instructions if necessary.
- › Completely empty the fluid container.
- › Thoroughly rinse the fluid container with water and dry.

## 12.4 Cleaning the sleeve part and adapter of the handpiece hose

- › Pull off the handpiece hose from the handpiece.
- › Unscrew and remove the sleeve part.

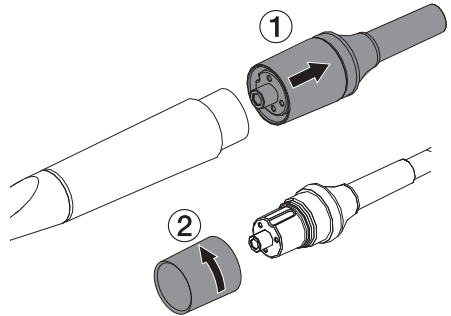

- › Clean the sleeve part and the adapter of the handpiece hose with a hygienic, soft brush and a moist, lint-free cloth.
- › Screw sleeve back in place.
- › Push the handpiece hose onto the handpiece.

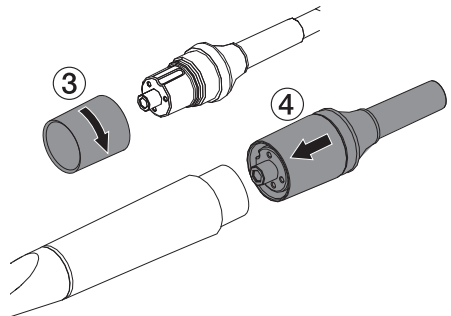

## 13 Reprocessing

### 13.1 Risk analysis and categorisation

A risk analysis and categorisation of medical products often used in dentistry must be performed before their reprocessing by the operator. Comply with all national directives, standards and specifications such as e. g. the "Recommendations from the Commission for Hospital Hygiene and Infection Prevention".

Accessories of the medical device are also subject to reprocessing.

Classification recommendation given intended use of the product: **semi-critical B to critical B**  
The operator is responsible for correct classification of the medical products, defining the reprocessing steps and performing the reprocessing.

### 13.2 Reprocessing procedure in accordance with EN ISO 17664

The reprocessing procedure after each patient treatment is carried out according to the reprocessing procedure established by EN ISO 17664.



#### Important information!







The reprocessing notes in accordance with EN ISO 17664 have been independently tested by Dürr Dental for the preparation of the device and its components for their reuse.

The person conducting the reprocessing is responsible for ensuring the reprocessing performed using the equipment, materials and personnel achieves the desired results. This requires validation and routine monitoring of the reprocessing process. Any deviation from the instructions described herein by the staff preparing the equipment could lead to lower effectiveness and possible negative consequences: these lie solely with the staff responsible.

Frequent reprocessing has little effect on the device components. The end of the product life cycle is especially influenced by the amount of wear and tear or damage resulting from its use.

The use of soiled, contaminated and damaged components is at the sole responsibility of the person performing the reprocessing and the operator.

The reprocessing method was validated as follows:

- **Pre-cleaning**  
  - FD 350 disinfection wipes (Dürr Dental)
  - Cleaning brush
- **Manual cleaning**  
  - ID 215 enzymatic instrument cleaner (Dürr Dental)
  - Cleaning brush
- **Manual disinfection**  
  - ID 213 instrument disinfection (Dürr Dental)
- **Automatic cleaning and disinfection**  
Was performed in accordance with EN ISO 15883 with tested efficacy.
  - Cleaning agent: Neodisher MediClean Forte
  - Washer-disinfector: PG 8535 (Miele)
  - Programmes: "Cleaning without neutralisation" and "THERMAL DES"
  - Rinsing adapter: Miele 68551101 D
  - Cleaning brush
- **Steam sterilisation**  
was performed in accordance with EN ISO 17665 with the fractionated vacuum procedure.
  - Pre-vacuum: 3 x
  - Sterilisation temperature: 132 °C
  - Sterilisation time: 4 minutes
  - Drying time: min. 20 minutes
- **Cleaning brush**  
Cleaning brush with nylon hairs, double-sided
  - Number of brush heads: 2
  - Brush material: nylon
  - Brush head length: 25 and 35 mm
  - Brush length: 5 and 10 mm

Example: Interlock cleaning brush, double-sided, green REF 09098

- › Comply with all national directives, standards and specifications for the cleaning, disinfection and sterilisation of medical products as well as the specific specifications for dental practices and clinics.
- › Comply with the specifications (see "13.6 Manual cleaning, intermediate rinsing, disinfection, final rinse, drying in ultrasonic bath" and "13.7 Automatic cleaning, intermediate rinsing, disinfection, final rinse, drying") when selecting the cleaning and disinfectant agents to be used.
- › Comply with the concentration, temperature, residence time and post-rinsing specifications issued by the manufacturer of the cleaning and disinfectant agent.
- › Only use cleaning agents that are non-fixing and aldehyde-free and display material compatibility with the product.
- › Only use disinfectants that are aldehyde-free and display material compatibility with the product.
- › Do not use any rinse aid (danger of toxic residue on the components).
- › Only use freshly-produced solutions.
- › Only use distilled or deionised water with a low bacterial count (at least drinking water quality) that is free from facultatively pathogenic microorganisms (e.g. legionella bacteria).
- › Use clean, dry, oil and particle-free compressed air.
- › Do not exceed temperatures of 138 °C.
- › Subject all devices used (ultrasonic bath, cleaning and disinfection device (CD), sealing device, steam steriliser) to regular maintenance and inspections.

## General information




### NOTICE


#### Equipment damage due to unsuitable products


Oils and care products containing oil will damage the device.

- › The handpiece must not be maintained with oil or with a care system that contains oil.

## EN 13.3 Preparation at the operating location

 Wear protective gloves.

 Wear protective goggles.

 Use a mask.

 Use protective clothing.



### WARNING

#### Risk of infection from contaminated products


Danger of cross contamination

- › Reprocess the product correctly and promptly before its first use and after every subsequent use.


### Rinsing the handpiece with water

Rinsing with water will rinse any remaining Vector polish residue from the handpiece and prevent blockages.

- › Fill the fluid container to approx. 1/3 with water.

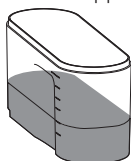


- › Place the handpiece in a sink or in a suitable container.

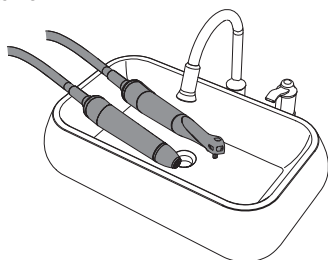

- › Start the rinsing process:  
Press the RINSE button for at least 2 seconds.

### RINSE

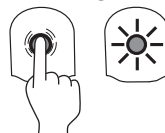


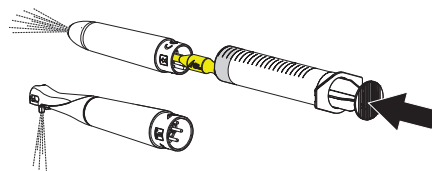
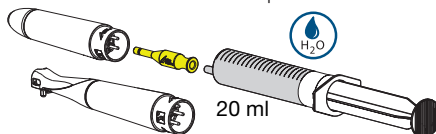
- › The LED flashes – the rinsing process takes around 30 seconds and ends automatically.

### Pre-cleaning



Do pre-cleaning of handpiece and accessories no more than 15 minutes after the unit has been used.

- › Clean the exterior surfaces completely with two cleaning cloths . Make sure that the surfaces are sufficiently moistened.
- › Note the action time of the cleaning agent.
- › Perform the procedure twice.
- › Draw 3 x 20 ml cold water (temperature < 20 °C) into a conventional sterile 20-ml disposable pipette with Luer connection and rinse the inner lumen of the handpiece.




### Transport

- › Protect the device from contamination when you transport it from the treatment chair to the reprocessing location.



### 13.4 Dismantling the handpiece

- › Unscrew the instrument, see "Inserting/changing instruments".
- › Take off the removable parts of the handpiece, Paro handpiece see "Disassembly", Scaler handpiece see "Disassembly".  
Unscrew union nut.

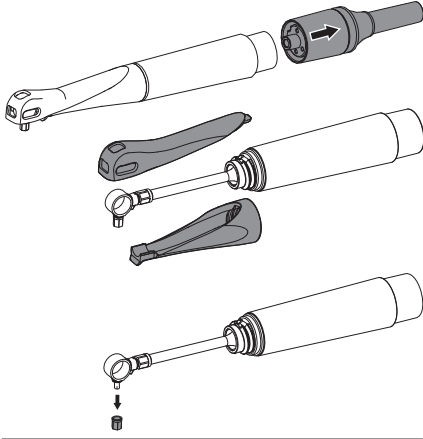

- › Check the Scaler handpiece light conductor for its light transmission, replace if necessary.

### 13.5 Manual cleaning, intermediate rinsing, disinfection, final rinsing, drying in the cleaning and disinfection bath


A disinfectant or combined cleaning and disinfectant agent is required for manual disinfection. It must have the following properties:

- certified, possibly virucidal efficacy (DWW/RKI, VAH or European Standards)


For further information, see: "General information".

### Cleaning



- › Place the removable parts of the handpiece (ring cover, rotary adaptor, scaler cover, light conductor), instrument holders of the tool kits (without instruments), torque wrench and the disassembled handpieces in the cleaning bath for the required action time, so that all parts are covered.

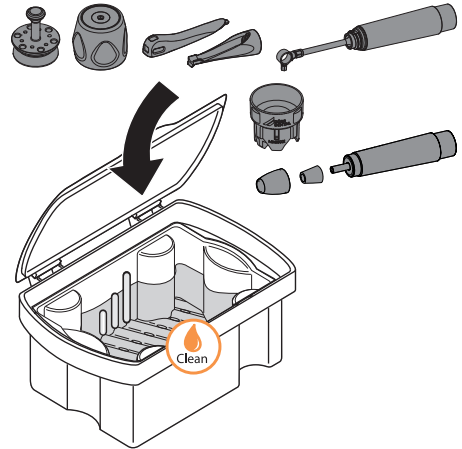



- › Clean all accessible external surfaces and internal surfaces below the surface for 5 minutes with a hygienic cleaning brush until all visible soiling has been removed.
- › Place instruments in small parts baskets into the cleaning bath.

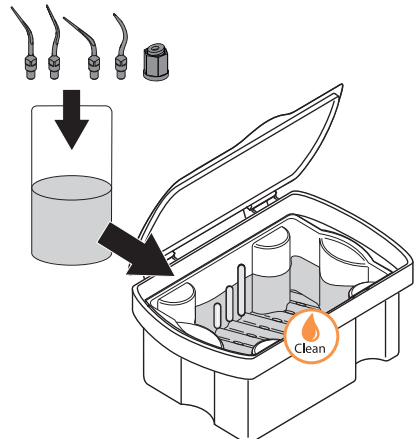

- › Rinse handpieces through at least 3 times using a 20-ml disposable pipette.

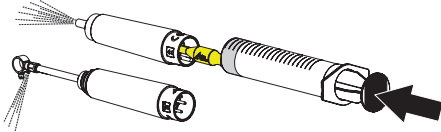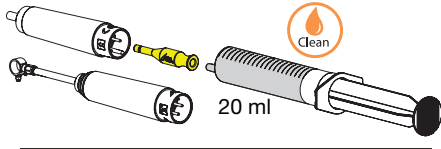

- › Screw Scaler instruments onto the rinsing adapter one after another and rinse through each of the instruments at least 3 times using a 20-ml disposable syringe.

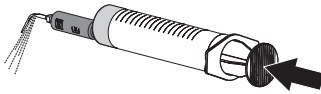
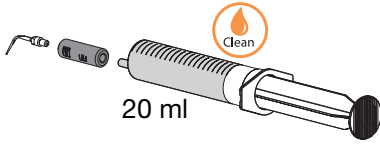


- › Remove all rinse adapters.
- › Comply with the cleaning agents' action times specified by the manufacturer.



### Intermediate rinsing

After the action time prescribed by the manufacturer:


- › Rinse all components under running water for at least 1 minute (temperature < 20 °C).
- › Rinse handpieces through with water at least 3 times using a 20-ml disposable pipette.

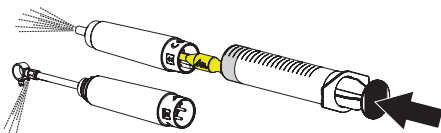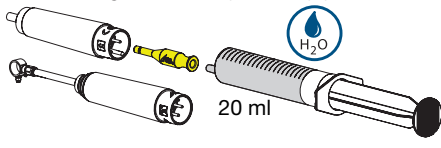

- › Screw scaler instruments onto the rinsing adapter one after another and rinse through each of the instruments at least 3 times with water using a 20-ml disposable pipette.

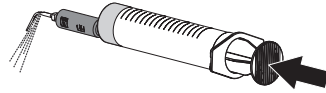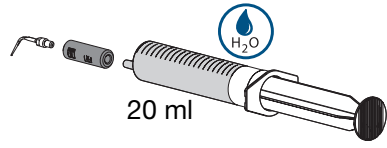


- › Areas that are difficult to reach, e.g. the Paro handpiece instrument chuck, should be thoroughly rinsed (at least 5 x for 5 seconds each) using an air and water syringe.

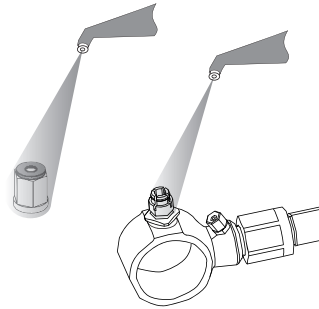


## Disinfection



- › Place the removable parts of the handpiece (ring cover, rotary adaptor, scaler cover, light conductor), instrument holders of the tool kits (without instruments), torque wrench and the disassembled handpieces in the disinfection bath for the required action time, so that all parts are covered.

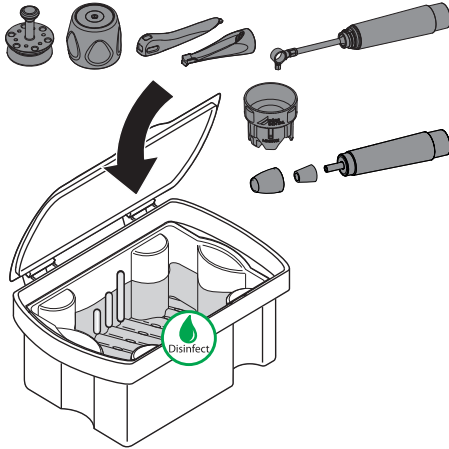



- › Place instruments in small parts baskets into the disinfection bath.

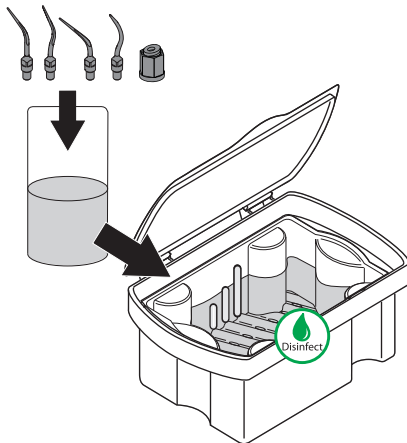

- › Rinse handpieces through at least 3 times using a 20-ml disposable pipette.

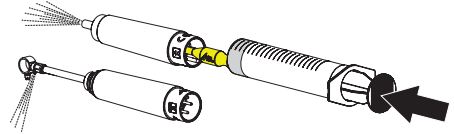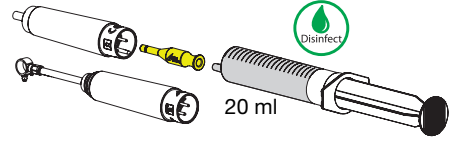

- › Screw Scaler instruments onto the rinsing adapter one after another and rinse through each of the instruments at least 3 times using a 20-ml disposable syringe.

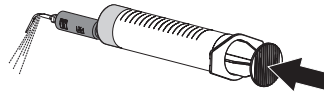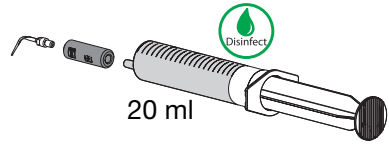

- › Remove all rinse adapters.
- › Follow the action times stated by the manufacturer for the cleaning agents and disinfectants.

## Final rinse

After the action time prescribed by the manufacturer:

- › Rinse all components under running water for at least 1 minute (temperature < 20 °C).

## Drying

- › If necessary, re-dry at a clean location using a hygienic, lint-free cloth.
- › Blow dry the components with compressed air in a clean location.

### 13.6 Manual cleaning, intermediate rinsing, disinfection, final rinse, drying in ultrasonic bath

A combined cleaning and disinfectant agent is required for manual cleaning and disinfection. It must have the following properties:

- certified, possibly virucidal efficacy (DW/RKI, VAH or European Standards)
- without chlorine, solvents, strong alkaline solutions (pH >11) or oxidising agents

For further information, see: "General information".

#### Cleaning in an ultrasonic bath




#### NOTICE

**Malfunctions in the handpiece due to improper handling during cleaning or disinfecting**


- > The handpieces of the Vector must only be cleaned or disinfected in a suitable container in an ultrasonic bath.
- > The handpieces must NOT be fully submerged in fluid.

- > Place the removable parts of the handpiece (ring cover, rotary adaptor, scaler cover, light conductor), instrument holders of the tool kits (without instruments) and torque wrench in the ultrasonic bath for the required action time, so that all parts are covered.

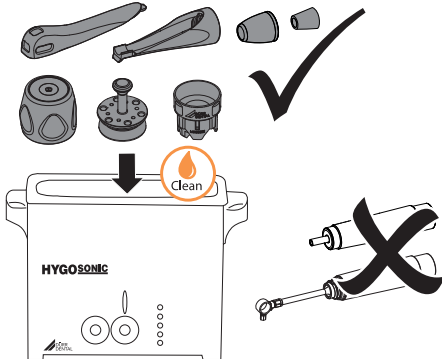


- > Clean all accessible external surfaces and internal surfaces below the surface for 5 minutes with a hygienic cleaning brush until all visible soiling has been removed.

- > Place instruments in small parts baskets into the ultrasonic bath.


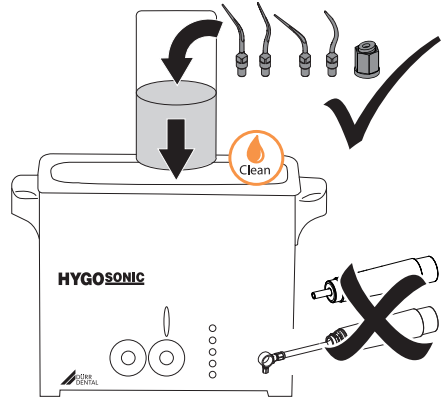


- > Place the handpiece without covers in a container with fluid. The drive mechanism of the handpiece must not lie in the fluid (malfunction). Therefore, observe the maximum fluid level for Paro and Scaler handpieces.

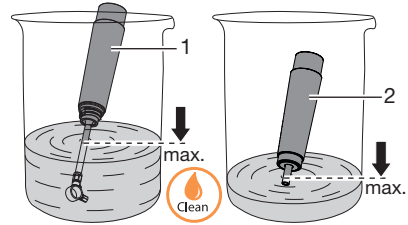



- 1 Paro handpiece
- 2 Scaler handpiece

- > Place container with the handpieces in the ultrasonic bath with a suitable carrier.
- > Rinse handpieces through at least 3 times using a 20-ml disposable pipette.

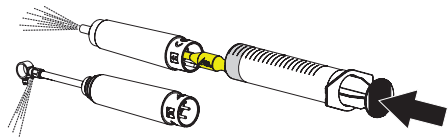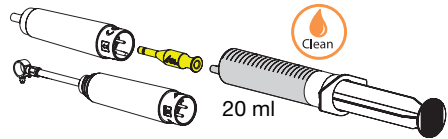

- > Remove all rinse adapters.

- › Screw Scaler instruments onto the rinsing adapter one after another and rinse through each of the instruments at least 3 times using a 20-ml disposable syringe.

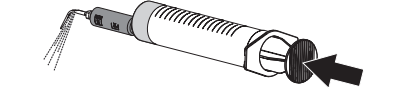
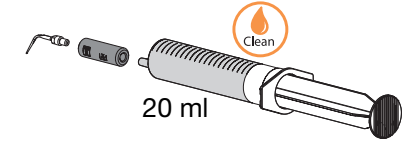


- › Comply with the cleaning agents' action times specified by the manufacturer.



### Intermediate rinsing

After the action time prescribed by the manufacturer:


- › Rinse all components under running water for at least 1 minute (temperature < 20 °C).
- › Rinse handpieces through with water at least 3 times using a 20-ml disposable pipette.

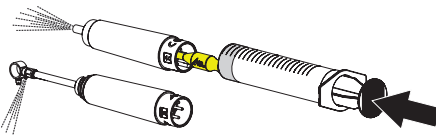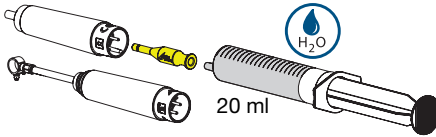

- › Screw scaler instruments onto the rinsing adapter one after another and rinse through each of the instruments at least 3 times with water using a 20-ml disposable pipette.

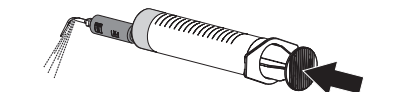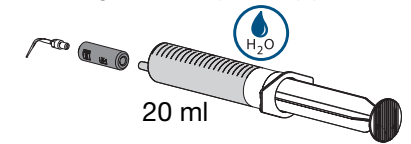

- › Areas that are difficult to reach, e.g. the Paro handpiece instrument chuck, should be thoroughly rinsed (at least 5 x for 5 seconds each) using an air and water syringe.

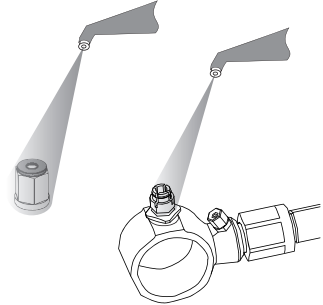

### Disinfection in an ultrasonic bath




#### NOTICE


#### Malfunctions in the handpiece due to improper handling during cleaning or disinfecting

- › The handpieces of the Vector must only be cleaned or disinfected in a suitable container in an ultrasonic bath.
- › The handpieces must NOT be fully submerged in fluid.


- › Place the removable parts of the handpiece (ring cover, rotary adaptor, scaler cover, light conductor), instrument holders of the tool kits (without instruments) and torque wrench in the ultrasonic bath for the required action time, so that all parts are covered.

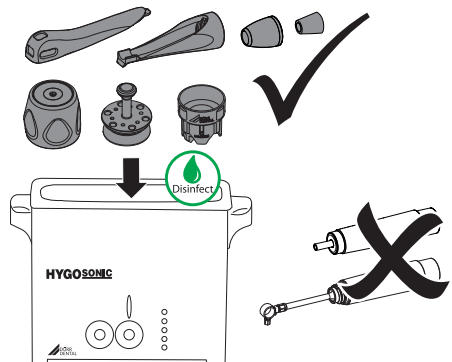

- › Place instruments in small parts baskets into the ultrasonic bath.


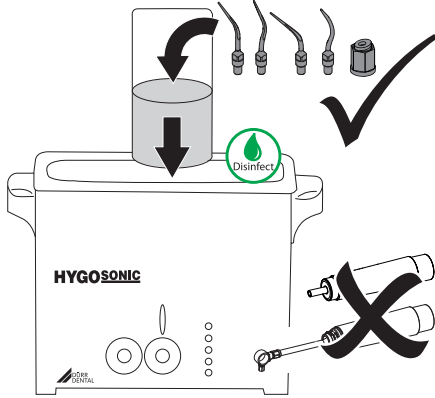


- › Place the handpiece without covers in a container with fluid. The drive mechanism of the handpiece must not lie in the fluid (malfunction). Therefore, observe the maximum fluid level for Paro and Scaler handpieces.


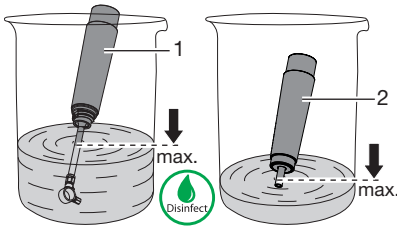


- 1 Paro handpiece
- 2 Scaler handpiece

- › Place container with the handpieces in the ultrasonic bath with a suitable carrier.
- › Rinse handpieces through at least 3 times using a 20-ml disposable pipette.

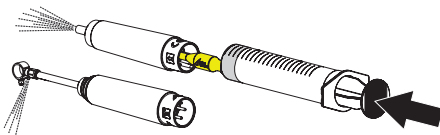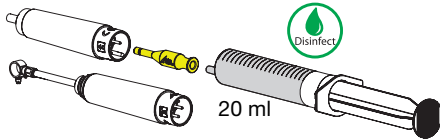

- › Screw Scaler instruments onto the rinsing adapter one after another and rinse through each of the instruments at least 3 times using a 20-ml disposable syringe.

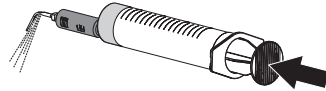
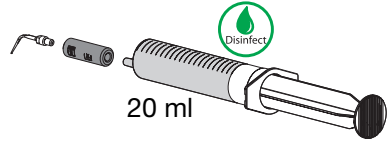


- › Remove all rinse adapters.
- › Follow the action times stated by the manufacturer for the cleaning agents and disinfectants.



### Final rinse

After the action time prescribed by the manufacturer:


- › Rinse all components under running water for at least 1 minute (temperature < 20 °C).
- › Rinse handpieces through with water at least 3 times using a 20-ml disposable pipette.

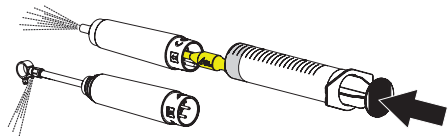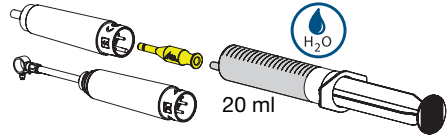

- › Screw scaler instruments onto the rinsing adapter one after another and rinse through each of the instruments at least 3 times with water using a 20-ml disposable pipette.

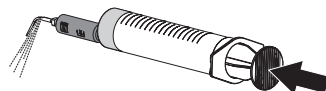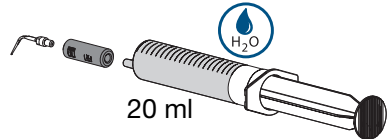

- › Areas that are difficult to reach, e.g. the Paro handpiece instrument chuck, should be thoroughly rinsed (at least 5 x for 5 seconds each) using an air and water syringe.

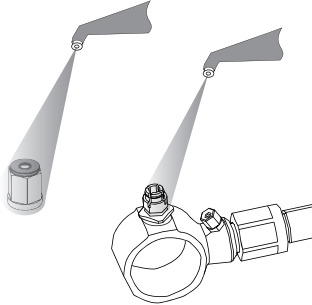

### Drying

- › If necessary, re-dry at a clean location using a hygienic, lint-free cloth.
- › Blow dry the components with compressed air in a clean location.

## 13.7 Automatic cleaning, intermediate rinsing, disinfection, final rinse, drying

### Selection of the washer-disinfector

Automatic cleaning and disinfection requires a washer-disinfector with the following properties and validated processes:

- Corresponds to and tested in accordance with EN ISO 15883
- Certified program for thermal disinfection ( $A_0$  value  $\geq 3000$  or at least 5 minutes at  $93^\circ\text{C}$ )
- Programme is suitable for the components and provides sufficient rinsing cycles.  
For more information: "General information".

### Selection of the machine cleaning agents and disinfectants

The following properties are required:

- Material compatibility with the product
- Corresponds with the manufacturer's specifications of the CD


For further information, see: "General information".

### Automatic cleaning and disinfecting



When arranging the parts in the washer-disinfector, make sure there are no areas missed by rinsing.

- › Place the handpiece on the special mountings for transmission instruments (e. g. Miele: ADS 2,  $\varnothing$  approx. 16 mm) in the washer-disinfector.

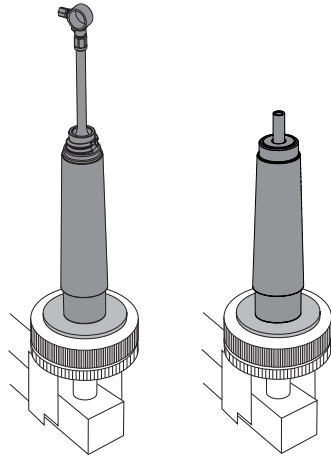

- › Place scaler instruments onto the special mountings for instruments (e.g. Miele: A 814) in the washer-disinfector.
- › Place Paro instruments in the instrument tray and place in the basket for small parts.
- › Fix removable parts of the handpiece (ring cover, rotary adaptor, scaler cover, light conductor), instrument holders of the tool kits (without instruments) and torque wrench with a suitable washer-disinfector holder.

## 13.8 Check for function

- › After the end of the cleaning and disinfection cycle, check the components for any residual soiling and moisture. If necessary, repeat the cycle.
- › Check the components for damage and replace if necessary.
- › The parts should be packaged as soon as possible after drying and checking.

## 13.9 Packing



### CAUTION

#### Endangering the sterilisation success

The fitted components are not reached by the steam and as such are not sterilised.

- › Do not fit the components before packaging.

For packaging of the components, use only sterile barrier systems made of transparent paper film that are approved for use in steam sterilisation according to the manufacturer information. This includes:

- Temperature resistance up to 138°C
- Standards DIN EN ISO 11607-1/2
- The applicable sections of the standard series DIN EN 868

The sterile barrier system must be large enough. Once it is loaded, the sterile barrier system must not be under any strain.

## 13.10 Steam sterilising



### WARNING

#### Health risk due to incorrect sterilisation

If the sterilisation not performed correctly, it may not be effective. The use of instruments that have not been properly sterilised can pose a health risk to the patient.

- › Only steam sterilisation must be used.
- › Comply with all of the specified process parameters.
- › Comply with the manufacturer's instructions regarding use of the steam steriliser.
- › Do not use any other methods.



### NOTICE

#### Damage to equipment due to incorrect sterilisation

If the sterilisation process is not performed correctly, this can cause damage to the product.

- › Comply with the manufacturer's instructions regarding use of the steam steriliser.
- › Comply with all of the specified process parameters.

#### Requirements placed on the steam steriliser:

- Corresponds to EN 13060 or EN 285 and/or ANSI AAMI ST79
- Suitable programme for the products listed (e. g. with hollow bodies, fractionated vacuum procedure in three vacuum steps)
- Sufficient product drying
- Validated process in accordance with DIN EN ISO 17665 (valid IQ/OQ and product-specific performance appraisal (PQ))

Perform the following steps:

- › Sterilise the parts for sterilisation (at least 20 minutes at 121°C, at least 4 minutes at 132°C or at least 5 minutes at 134°C).



Do not exceed 138 °C.

#### Marking

- › Mark the packaged, treated medical product in such a way as to ensure safe application.

## 13.11 Issue clearance for the parts for sterilisation

The reprocessing of the medical products ends with the documented clearance for storage and renewed use.

- › Document the clearance of the medical product after reprocessing.



## 13.12 Storing parts for sterilisation

- › Comply with the stated storage conditions:
  - Store the parts protected against contamination
  - Dust-protected, e.g. in a locked cabinet
  - Protected against moisture
  - Protected against excessive temperature fluctuations
  - Protected against damage

Packaging for a sterile medical device can suffer damage as a result of a particular incident and the passage of time.

Potential external contamination of the sterile barrier system should be taken into account in terms of aseptic preparation when establishing the storage conditions.

## 14 Treatment breaks for more than 24 hours

If no treatment is carried out for a period of 24 hours or more, reprocessing of the hose system must be performed.


### 14.1 Cleaning and disinfecting the hose system

The hose system is disinfected using the ready-to-use, aldehyde-free solution *Vector/RinsEndo Disinfectant*.

**Preparation:**


**i** Before starting the CLEAN or RINSE (rinsing/disinfecting) functions, if there is a fluid bag in the unit it will need to be removed. If an attempt is made to start the function with a fluid bag inserted, a long, low audible signal will sound and the action will not be performed.

- › Remove the cover of the fluid bag.
- › Push the piercing mechanism to the rear.
- › Take out the fluid bag.
- › Slide the piercing mechanism forwards as far as it will go.

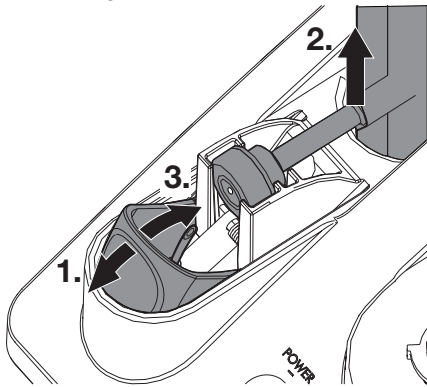


- › Pull off the hose connection from the handpiece.

- › Lay the handpiece hose in the sink.

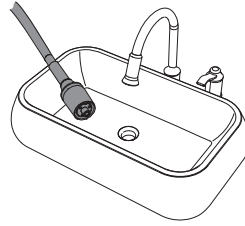


**Rinsing with water:**

- › Fill the fluid container to approx. 1/3 with water.

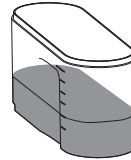

- › Touch the RINSE (rinsing/disinfecting) button for at least 2 seconds.


**RINSE**

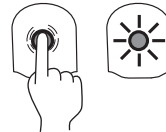


- › LED flashes – the rinsing process takes around 30 seconds and ends automatically.
- › Rinsing with water will flush out any remaining *Vector Polish* residue.
- › After the end of the rinsing process drain any remaining fluid in the system.

**Disinfection with *Vector/RinsEndo Disinfectant*:**

- › Pour 2 sealing caps (approx. 40 ml) of the *Vector/RinsEndo Disinfectant* in the fluid container.

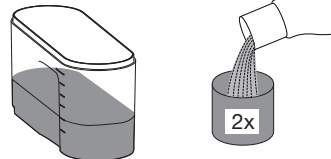

- › Touch the RINSE (rinsing/disinfecting) button for at least 2 seconds.

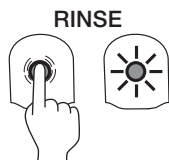

- › LED flashes – the *Vector/RinsEndo Disinfectant* is rinsed into the system and the process ends automatically. The *Vector/RinsEndo Disinfectant* stays in the system until the next treatment.
- › Empty any remaining *Vector/RinsEndo Disinfectant* from the fluid container.
- › Thoroughly rinse the fluid container with water and dry.

- › Fluid pre-delivery to the handpiece: Touch the FLUID button for at least 2 seconds. The LED flashes simultaneously with the LED in the piercing mechanism; the process ends automatically.

### 2. *Reprocessing not performed before the treatment break:*

- › Perform reprocessing before initial start-up of the unit, "13 Reprocessing".

## 14.2 Initial start-up after a break in treatment for more than 24 hours

The initial start-up process will depend upon whether the hose system was reprocessed before the break in treatment. Proceed as follows depending on the situation:

### 1. *Reprocessing performed before the treatment break:*

- › Rinse the system with water.



Thorough rinsing with water will flush out any disinfectant remaining in the hoses, and any irritation to the patient caused by the taste of residual disinfectant will be avoided.

- › Pull off the hose connection from the handpiece.
- › Lay the handpiece hose in the sink.
- › Fill the fluid container with water.
- › Touch the RINSE (rinsing/disinfecting) button for at least 2 seconds.  
LED flashes – the rinsing process takes around 30 seconds and ends automatically.

### *Preparing the device for treatment:*

- › Fill the fluid container.
- › Insert a fluid bag.

## 15 Maintenance

### 15.1 Service kit

The service kit has been configured for the Vector System.

The service kit contains spare parts and tools, which can be used for maintenance and repairs.

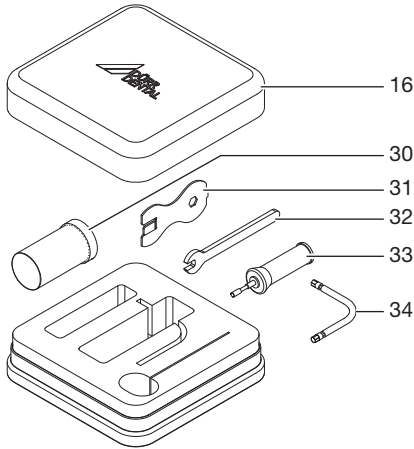


Fig. 1: Service kit

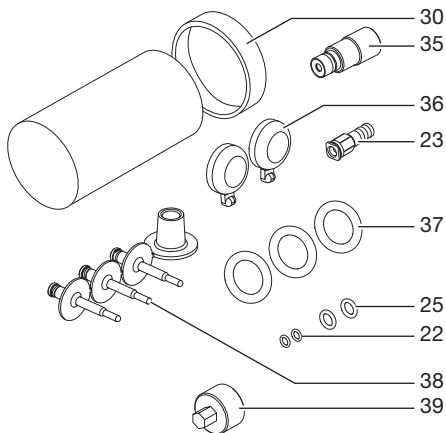



Fig. 2: Container for small parts with content

- 16 Service kit
- 22 O-rings
- 23 Interchangeable bushing, mounted
- 25 Seals, green
- 30 Container for small parts
- 31 Combination wrench for sealing screw

- 32 Flat wrench SW5
- 33 Nozzle torque wrench
- 34 Test tool, SW3.5, for torque wrench in the tool kit cover
- 35 Sealing screw
- 36 Rubber seal in piercing mechanism
- 37 O-Rings for valve in fluid container\*
- 38 Assembly device for the support seal of the nozzle
- 39 Adapter for tool kit

### 15.2 Changing the seal


- › Change the rubber seal in the piercing mechanism when visible signs of wear become apparent.
- › Replace the rubber seal in the piercing mechanism immediately if it is missing.

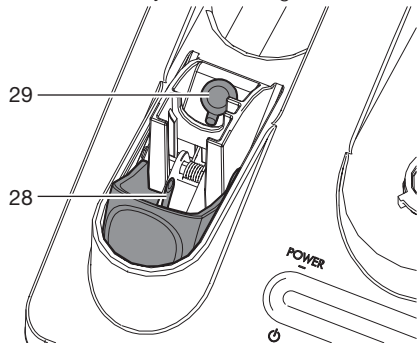

- 28 Piercing mechanism
- 29 Rubber seal

### 15.3 Changing the valve in the fluid container


The valve on the underside of the fluid container must be cleaned regularly and checked for blockages and leaks.

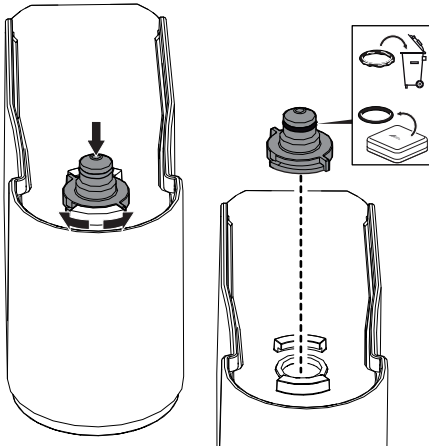

- › Press the valve lightly against the fluid container and unscrew it anti-clockwise.
- › Clean the valve.  
If cleaning is not possible, e.g. if the filter in the valve is blocked, then the valve needs to be replaced.
- › Check the O-ring.  
In the event of leaks, loose seating or visible damage the O-ring must be replaced.
- › Insert the valve in the holder and tighten it clockwise as far as it will go.

### 15.4 Checking instrument wear

Instrument wear can be checked using the Vector toolcard:

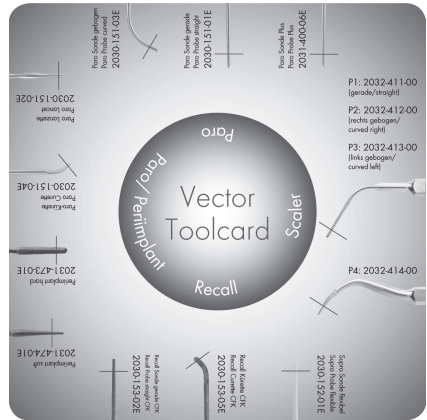

- Place the handpiece against the toolcard. If the instrument tip extends beyond the red marking then the instrument can still be used.
- If the instrument tip just reaches the red marking then the instrument shows signs of wear but can still be used.
- If the instrument tip does not reach the red marking then the instrument must be replaced.


## 15.5 Replacing the light conductor in the Scaler handpiece

The light conductor must be checked regularly for light transmission. Over the course of time it becomes opaque or takes on a milky colour. This impairs its ability to function and means that it needs to be replaced.



The light conductor can be sterilised several times. As soon as it starts to become opaque or take on a milky colour, light transmission begins to decrease.

- › Unscrew the cover.
- › Disconnect the light conductor.
- › Connect a new light conductor.
- › Screw on the cover.

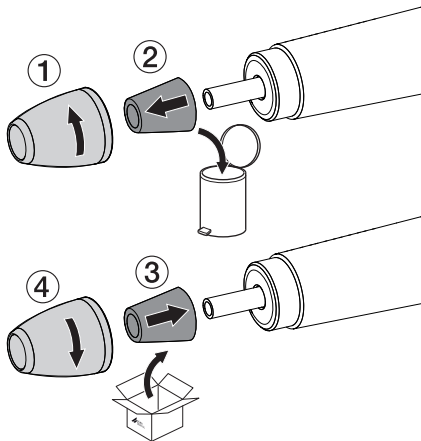

## 15.6 Replacing the sealing screw

The sealing screw in the adapter of the handpiece hose must be replaced 1x per year.

### Unscrew the sealing screw:


- › Unscrew the sealing screw by hand. If the screw is stuck, use the combination wrench from the service kit.

### Screw in the sealing screw:



Never use the combination wrench to screw in the sealing screw. If it is screwed in with too much torque, this can cause damage to the mixing chamber.

- › Firmly tighten the sealing screw by hand.

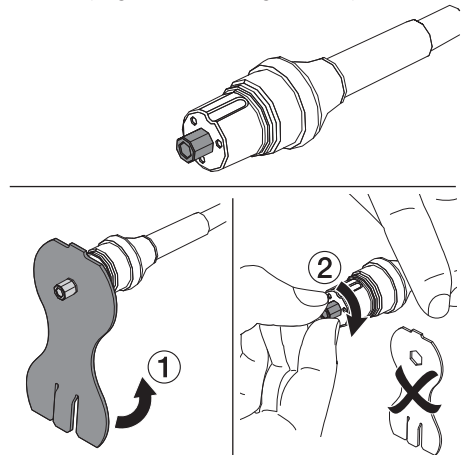


## 15.7 Inserting or changing the battery in the flexible foot switch

A new battery should be inserted in the flexible foot switch prior to initial start-up in wireless

operation or if the power of the existing battery is low.


### *Opening the cover:*

- › Press both pins on the flexible foot switch together at the same time and lift the cover off.

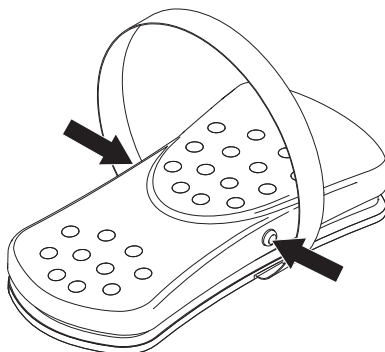

### *Checking the battery:*

- › Press the button to the left of the battery. Green LED lights up: battery power is OK. Green LED does not light up: change the battery.

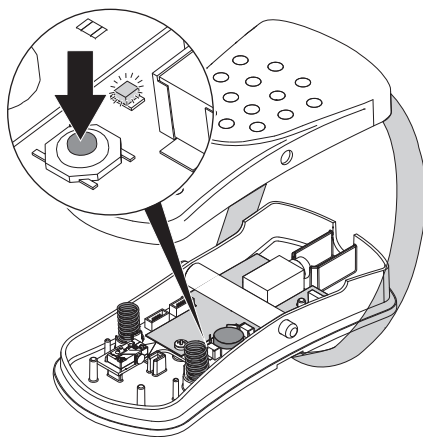


### *Taking out the battery:*




The old battery must be properly disposed of in accordance with applicable national and regional guidelines for environmentally friendly disposal. Do not throw batteries away with domestic waste.

- › Take out the battery from the holder and dispose of it in an environmentally friendly way.

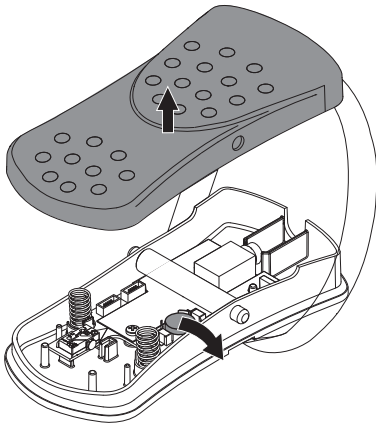

**Inserting the battery:**


- › Insert the battery in the holder. Make sure that the battery is inserted with the correct polarity.

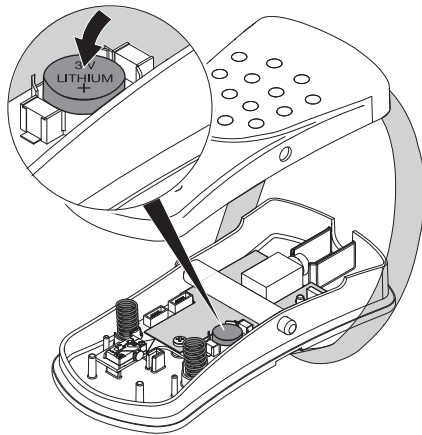

**Closing the cover:**

- › Make sure that the two pedal return springs are present and correctly positioned. If they are not, the function of the unit may be impaired.


- › Fit the cover in such a way that the two pins on the side of the flexible foot switch snap into the holes in the cover.

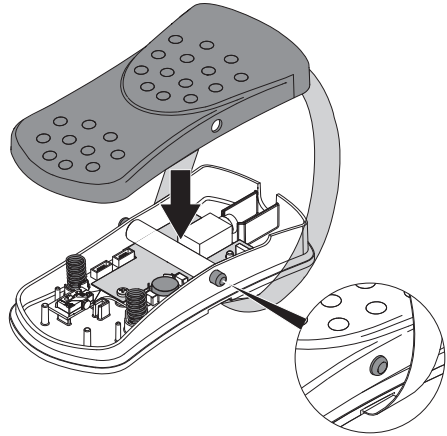




## 15.8 Changing the nozzle in the Paro handpiece

- › Remove the instrument, "Inserting/changing instruments".

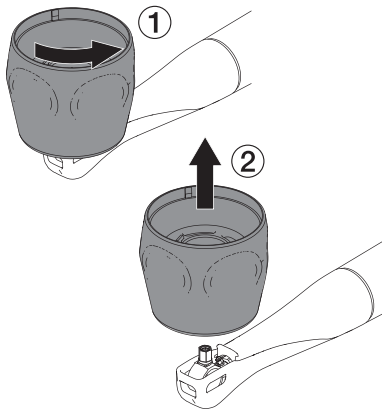

- › Place the nozzle torque wrench (contained in the service kit) on the nozzle and unscrew the nozzle anti-clockwise.

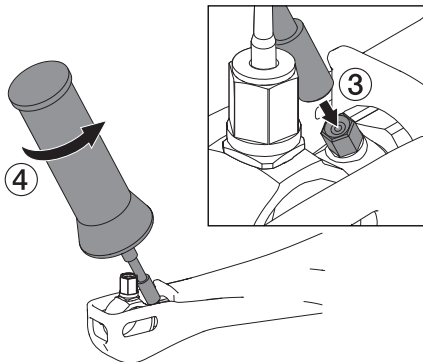



- › Rinse the thread in the resonating body with the air and water syringe.



New nozzles are stored in the nozzle torque wrench housing.

The nozzles are made of plastic and do not have a thread. The thread die in the resonating body is self-tapping.

- › Insert the new nozzle with the hexagon into the nozzle torque wrench and, while maintaining the correct angular position, apply light force to screw it into the thread cutting die of the resonating body.

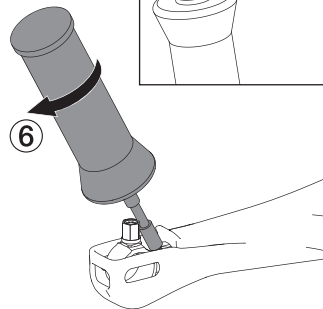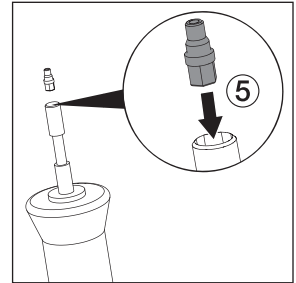


### WARNING

#### Risk of injury due to loose nozzles


Injuries can occur if the nozzle is not screwed in correctly and becomes loose during treatment.

- › Only work with the nozzle properly tightened. (An audible click can be heard when screwing in the nozzle).
- › As soon as the self-tapping thread takes hold, reduce the pressure. Screw on the nozzle until a slight cracking noise is heard (torque limit is reached). If it is not possible to tightly screw the nozzle into position (torque limit is not reached or the nozzle falls out again), check whether there is a chip in the thread. Remove the chip with a fine needle.

- › Check the nozzle function with an instrument clamped in place.


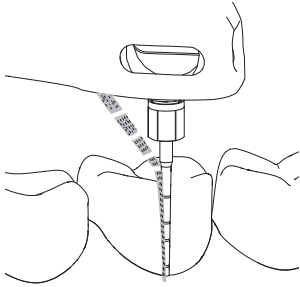


- › Carefully and completely remove the old or defective support seal from the union nut using a suitable tool.

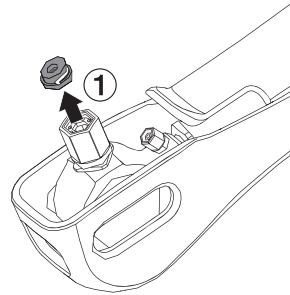

### 15.9 Changing the support seal


The support seal (supplied in the service kit) should be changed every 6 months or in the following situations:

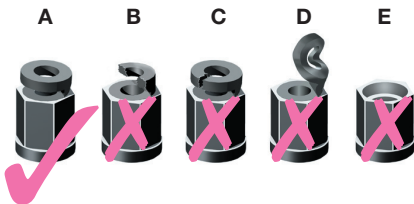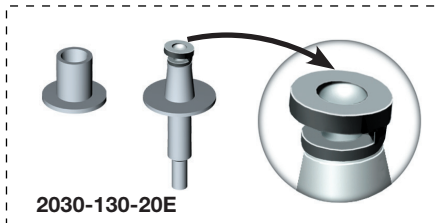


Condition of support seal:

- ✓ A – OK
- ✓ B – torn off
- ✓ C – cracked open
- ✓ D – twisted backwards
- ✓ E – missing

- › Place the assembly sleeve onto the union nut.
- › Insert the assembly pin with the sealing ring as shown in the picture.

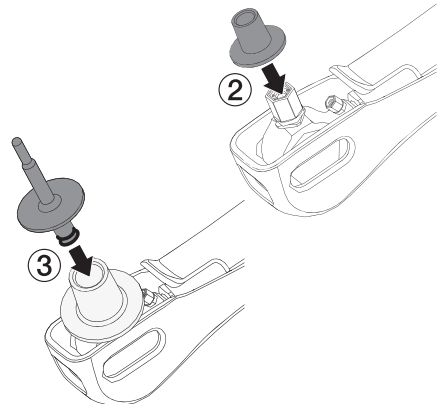

- › Press the assembly pin downwards until it is felt to engage.

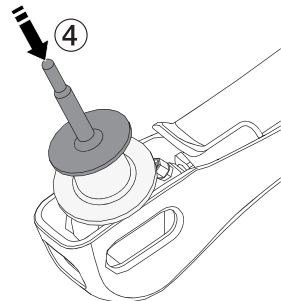

- › Remove the assembly sleeve and store it in a safe place for the next time it is needed.

**Result:**

The assembly pin is a disposable item and can be disposed of after use.



Disinfect the assembly sleeve with standard commercially available disinfectant, e.g. FD 322 or ID 212 forte. The assembly sleeve is not autoclavable.

### 15.10 Checking the function of the tool kit cover



The torque wrench located in the tool kit cover is subject to wear as a result of various factors. If the torque wrench is defective, this means that the tools in the handpiece can no longer be properly secured. For this reason the function of the torque wrench should be regularly checked.

- › Place the test tool (included in the service kit) vertically in the torque wrench.

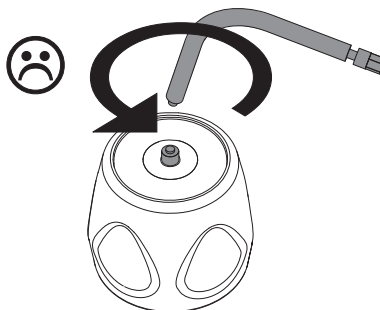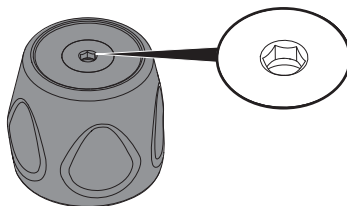

Never rotate the test tool anti-clockwise, as otherwise it will break off.

- › Rotate the test tool clockwise. Hold the tool kit cover steady while doing this.

**Result:**

Test tool stays intact:

- Tool kit cover can still be used further.

Test tool breaks off:

- Tool kit cover can no longer be used and must be replaced with a new one.

## 15.11 Changing the sealing ring of the resonant ring

Sealing ring, order number




### NOTICE

**Damage to the resonant ring due to twisting**


› When loosening and tightening the nut do not hold down the resonant ring.

- › Dismantle the handpiece, "Disassembly".
- › Place the sealing ring and SW 5 flat wrench from the service kit ready.

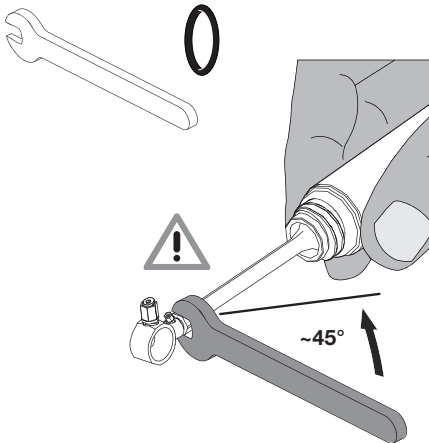

- › Loosen the nut using the flat wrench, turning approx. 45° (1/4 turn) anti-clockwise.

- › Change the sealing ring of the resonant ring.

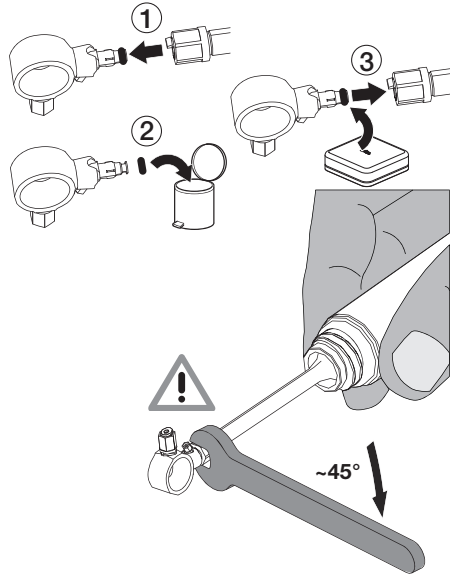

- › Tighten the nut using the flat wrench, turning approx. 45° (1/4 turn) clockwise.
- › Mount the handpiece, "Assembly".

## 15.12 Changing the resonant ring

Resonant ring, order number




### NOTICE

#### Damage to the resonant ring due to twisting

- › When loosening and tightening the nut do not hold down the resonant ring.

- › Dismantle the handpiece, "Disassembly".

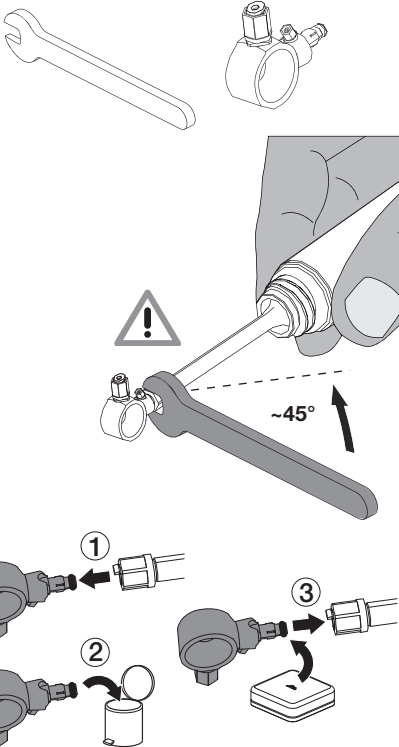

- › Get the SW5 flat wrench from the service kit and the resonant ring ready.
- › Loosen the nut using the flat wrench, turning approx. 45° (1/4 turn) anti-clockwise.
- › Replace the resonant ring.
- › Tighten the nut using the flat wrench, turning approx. 45° (1/4 turn) clockwise.
- › Mount the handpiece, "Assembly".

## 15.13 Changing the interchangeable bushing with union nut


The interchangeable bushing (instrument chuck) and union nut are subject to wear due to usage. They need to be replaced in the following cases:

- If there is visible damage.
- If the instrument cannot be securely attached.




The interchangeable bushing, union nut and mounting pin are contained in the service kit as an assembled unit.

- › Unscrew the union nut using the tool kit cover.

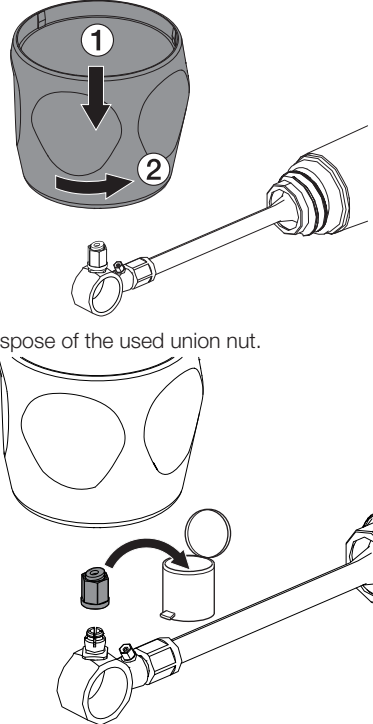


- › Dispose of the used union nut.


- › Place the adapter (contained in the service kit) in the hexagon socket of the tool kit cover and unscrew the interchangeable bushing from the resonant ring.

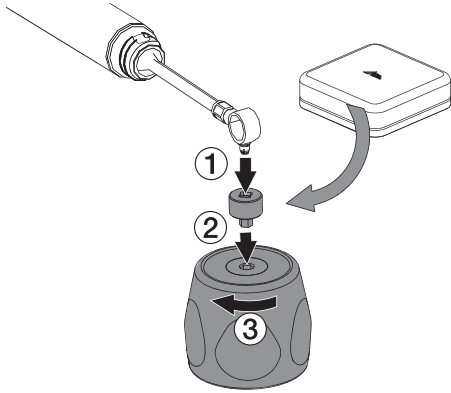


- › Dispose of the used interchangeable bushing.

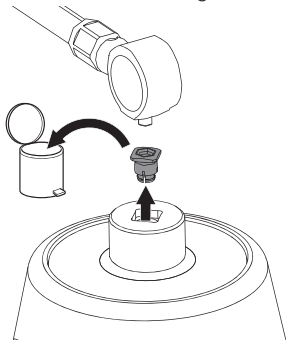



- › Tighten the interchangeable bushing with union nut and the mounting pin by hand on the resonant ring.

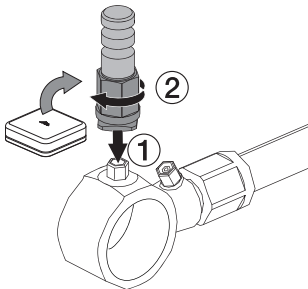


- › Place the adapter in the hexagon socket of the tool kit cover and tightly screw the interchangeable bushing onto the resonant ring until the torque limit is reached and ratcheting of the tool kit is clearly heard.

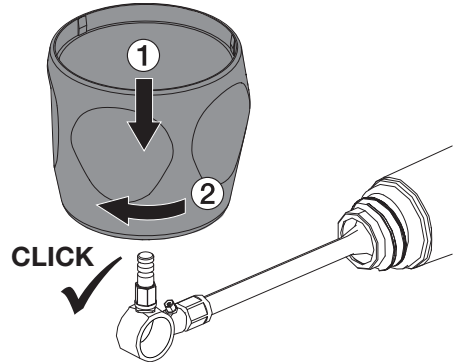

- › Loosen the union nut by turning it anti-clockwise through 90°. Remove and dispose of the assembly pin.

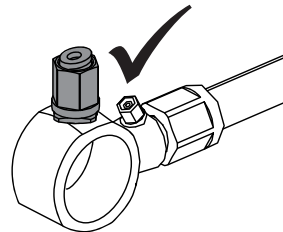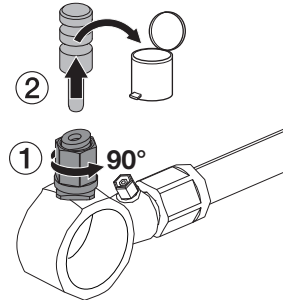

- › Take out the adapter from the tool kit cover and store it in the service kit.

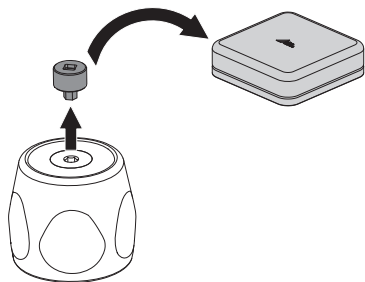

# ? Troubleshooting

## 16 Tips for operators and service technicians



Prior to working on the device or in case of danger, disconnect it from the mains (e. g. pull the mains plug).





Any repairs exceeding routine maintenance may only be carried out by qualified personnel or our service.

Error	Possible cause	Remedy
<b>Unit does not start</b>	Unit is not switched on.	› Touch the On / Standby button for at least 2 seconds.
	The flexible foot switch cable is not connected.	› Connect the flexible foot switch cable.
	The flexible foot switch is defective.	› Replace the flexible foot switch and send in the defective flexible foot switch for repairs.
	Handpiece is defective.	› Replace the handpiece and send in the defective handpiece for repairs.
	The resonating body in the Paro handpiece is defective.	› Refer to Troubleshooting, point 12.
<b>Unit runs with interruptions in wireless operation or wireless operation with the flexible foot switch is not possible.</b>	The battery in the flexible foot switch is empty.	› Check battery performance and insert a new battery if necessary.
	Interference caused by external wireless signals.	› Operate the flexible foot switch using the foot switch cable.
	The wireless module in the flexible foot switch is defective.	› Operate the flexible foot switch using the foot switch cable or call a Service Technician.
	The wireless module in the basic unit is defective.	› Operate the flexible foot switch using the foot switch cable or call a Service Technician.
	Pairing of the flexible foot switch has not been performed.	› Before using it for the first time, pairing (i.e. synchronisation/coupling) must be performed between the flexible foot switch and the unit.



Error	Possible cause	Remedy
<b>Operating the flexible foot switch does not activate the handpiece.</b>	Unit is not switched on.	› Switch on the unit.
	The connector of the flexible foot switch cable is not plugged in correctly.	› Plug in the connector correctly.
	The flexible foot switch cable is defective.	› Replace the foot switch cable.
<b>Fluid is sprayed in pulses onto the instrument (Vector Paro handpiece).</b>	Normal operating status. <b>NOT A FAULT.</b> The avoidance of spray mist and instrument warming means that only a minimal quantity of liquid is required for cooling.	
<b>Pulsed application of fluid is not clean, or fluid drips out afterwards</b>	Fluid bag has not been pierced, or the piercing mechanism is not fully closed.	› Pierce the fluid bag; to do this, fully close the piercing mechanism.
	Grey rubber seal is missing or defective.	› Mount a new rubber seal.
	Fluid container is empty.	› Fill the fluid container.
	The O-ring on the valve of the fluid container is not leak-tight.	› Replace the O-ring or the valve of the fluid container.
	Overall system is not correctly bled.	› Fill the fluid container with water. › Touch the RINSE (rinsing/disinfecting) button for at least 2 seconds.
<b>Fluid leaks between handpiece and handpiece hose during operation</b>	Handpiece is not correctly attached to the handpiece hose.	› Correctly attach the handpiece to the handpiece hose.
	Sealing screw in handpiece hose leaks.	› Pull off the handpiece from the handpiece hose.
		› Replacing the sealing screw.
<b>Fluid leaks out between Scaler handpiece and instrument.</b>	Scaler instrument was not screwed on with the full torque.	› Screw on the Scaler instrument correctly.
	Scaler instrument has become loose.	› Screw the Scaler instrument in place correctly.
		› If there are signs of wear, replace the instrument.
<b>Fluid leaks between handpiece and handpiece hose during operation.</b>	Air in the fluid system.	› Fill the container with water. › Touch the FLUID button for at least 2 seconds.
<b>Fluid leaks from the underside of the basic unit.</b>	Hose connections within the unit have become loose or are defective.	› Call a Service Technician.

Error	Possible cause	Remedy
<b>No fluid ejected from the nozzle of the Paro handpiece.</b>	The nozzle in the Paro handpiece is blocked or defective.	› Replace the nozzle in the Paro handpiece.
	Fluid container is empty.	› Fill the fluid container.
	Handpiece blocked.	› Clean the handpiece.
	Fluid pump is defective.	› Call a Service Technician.
<b>No fluid applied from the nozzle in the Paro handpiece when working with the flexible foot switch.</b>	Fluid release is switched off.	› Press the FLUID button.
	Fluid bag is empty.	› Change the fluid bag.
	No fluid bag inserted.	› Insert a fluid bag.
	Fluid bag has not been pierced, or the piercing mechanism is not fully closed.	› Pierce the fluid bag; to do this, fully close the piercing mechanism.
	Sealing screw on the handpiece hose is blocked or defective.	› Pull off the handpiece from the handpiece hose. › Replace the sealing screw (replacement part in the service kit).
<b>Fluid leaks at the coupling of the resonant body in the Paro handpiece.</b>	The sealing ring on the coupling of the resonating body is defective.	› Change the sealing ring of the resonant ring.
	The nozzle in the Paro handpiece is blocked or defective.	› Replace the nozzle in the Paro handpiece.
<b>Unusual noises from the Paro handpiece</b>	Union nut in the Paro handpiece is loose.	› Insert the instrument as far as it will go and tighten the union nut using the tool kit cover.
	Support seal in the union nut is missing or defective.	› Replace the support seal (replacement part in the service kit).
	Resonant ring is loose.	› Tighten the nut on the resonant ring.
	Resonant ring is defective.	› Replace the resonant ring.
<b>Paro handpiece ring cover does not engage correctly.</b>	Rotary adaptor not fully closed.	› Fully close the rotary adaptor as far as it will go.
	Ring cover is defective.	› Replace the ring cover.
<b>Handpiece cannot be attached to the handpiece hose.</b>	The O-ring of the sealing screw has dried out or is defective.	› Grease the O-ring and replace the sealing screw if required (replacement part in the service kit).
	Contact pins are bent.	› Send in the handpiece.
	Water connector is bent.	› Send in the handpiece.

Error	Possible cause	Remedy
<b>Instrument cannot be inserted or instrument is loose in the instrument chuck.</b>	Instrument chuck is deformed.  If the locking nut is tightened using the torque wrench without an instrument being in position, then the instrument chuck can be deformed as a result.	<ul style="list-style-type: none"> <li>› Unscrew the union nut using the tool kit cover.</li> <li>› Carefully widen the instrument chuck using an appropriate tool, e.g. "Supra probe flexible", until the instrument can be correctly inserted again.</li> </ul>
	Instrument chuck is dirty.	<ul style="list-style-type: none"> <li>› Unscrew the union nut using the tool kit cover.</li> <li>› Clean the union nut and the instrument chuck with the air and water syringe.</li> <li>› Check for signs of damage and for completeness.</li> </ul>
	Instrument chuck is worn.	<ul style="list-style-type: none"> <li>› Replace the interchangeable bushing with union nut.</li> <li> The instructions for replacement are enclosed with the set.</li> </ul>
<b>Scaler instrument cannot be loosened at all or only with difficulty.</b>	During working with polish, polish has collected in the thread between the Scaler instrument and the Scaler handpiece.	<ul style="list-style-type: none"> <li>› Place the Scaler handpiece with instrument in the active ultrasonic bath until it is covered and leave it standing for a few minutes.</li> <li>› Then loosen the Scaler instrument. If this is not possible, repeat the process.</li> </ul>
<b>Resonant ring has become twisted during operation (Paro handpiece)</b>	Loosening the union nut has broken the locking mechanism.	<ul style="list-style-type: none"> <li>› Replace the handpiece.</li> </ul>
	Handpiece is defective.	<ul style="list-style-type: none"> <li>› Replace the handpiece.</li> </ul>

Error	Possible cause	Remedy
<b>Orange "FUNCTION" LED lights up or flashes</b>	The treatment was interrupted as water may be present in the following areas: - between the handpiece and the handpiece hose. - at the instrument chuck of the Paro handpiece (union nut has loosened slightly).	<ul style="list-style-type: none"> <li>› Clean the affected areas and dry using the air and water syringe.</li> <li>› If necessary, tighten the union nut at the instrument chuck of the Paro handpiece.</li> </ul>
	Instrument contact pressure too high during treatment.	<ul style="list-style-type: none"> <li>› Reduce the contact pressure and press the flexible foot switch, the LED will go out.</li> <li>› If necessary, the instrument might require spraying and drying using the air and water syringe.</li> </ul>
	Interchangeable bushing or resonant ring of the Paro handpiece defective.	<ul style="list-style-type: none"> <li>› Replace the interchangeable bushing with union nut, order no.</li> <li>› Replace the resonant ring, order no. .</li> <li><b>i</b> The instructions for performing the replacement are enclosed with the sets.</li> </ul>
	Instrument defective.	<ul style="list-style-type: none"> <li>› Replace the instrument.</li> </ul>
	Handpiece defective	<ul style="list-style-type: none"> <li>› Replace the handpiece. Send in the defective handpiece for repairs.</li> </ul>
<b>Orange "PEDAL" LED lights up</b>	The power of the battery in the flexible foot switch is low.	<ul style="list-style-type: none"> <li>› Check battery performance and insert a new battery if necessary.</li> </ul>
<b>Orange "PEDAL" LED flashes,</b>	No flexible foot switch connected (cable operation) or paired (wireless operation).	<ul style="list-style-type: none"> <li>› Connect a flexible foot switch (cable operation) or perform pairing (wireless operation).</li> </ul>
<b>Blue "CLEAN" LED continues to light up after a cleaning program has been carried out.</b>	The CLEAN function of the cleaning process was not fully completed, or it was stopped early.	<ul style="list-style-type: none"> <li>› Carry out the cleaning process CLEAN in full.</li> </ul>
<b>The blue "POWER" and "LIQUID" LEDs fail to light up after the unit is switched on.</b>	Handpiece is not recognised by the basic unit.	<ul style="list-style-type: none"> <li>› Attach a different handpiece in position. If the "POWER" and "LIQUID" continue to fail to light up, call a Service Technician.</li> </ul>
	Handpiece is defective.	<ul style="list-style-type: none"> <li>› Attach a different handpiece in position. If the "POWER" and "LIQUID" continue to fail to light up, call a Service Technician.</li> </ul>

Error	Possible cause	Remedy
<p><b>The illumination in the Scaler handpiece becomes increasingly dimmer.</b></p>	<p>The light conductor has become opaque or has become milky.</p> <hr/> <p>Illumination LEDs are defective.</p>	<p>› Replace the light conductor.</p> <hr/> <p>› Send in the defective Scaler handpiece for repairs.</p>







**Hersteller/Manufacturer:**

DÜRR DENTAL SE  
Höfigheimer Str. 17  
74321 Bietigheim-Bissingen  
Germany  
Fon: +49 7142 705-0  
[www.duerrdental.com](http://www.duerrdental.com)  
[info@duerrdental.com](mailto:info@duerrdental.com)

