



**33** CLINICAL EVALUATORS

---

**657** TOTAL USES

---

**94%** CLINICAL RATING

## Description

**Aura Bulk Fill** is a nanohybrid composite designed for filling posterior teeth in a single layer. The Ultra High Density (UHD) glass filler provides a high-strength interface that can withstand high compressive forces while still delivering excellent handling and polishing properties.



## Unique Attributes

**Aura Bulk Fill** offers a 5 mm depth of cure and one universal shade to simplify posterior composite placement.

## Clinical Evaluation Results

- **Aura Bulk Fill** was evaluated by 33 consultants in 657 uses.
- The composite handles well, adapting to cavity walls and allowing sculpting without sticking to the instrument.
- The single, universal shade approximates A2 and blends well with surrounding enamel.
- Options of two extra shades were requested by some consultants.



## Results from THE DENTAL ADVISOR Biomaterials Research Laboratory

- Gloss after polishing ..... 83 (2) gloss units
- Radiopacity 2.9 (0.1) mm of aluminum per mm of material
- Depth of cure (20 s) ..... 7.3 mm
- Depth of cure (40 s) ..... 8.9 mm

## Evaluators' Comments


"I can finish most of my posterior restorations in a single layer because of the deep curing."

"Easy to extrude and handling is amazing."


"The universal shade is really just that. I did not expect the color to blend with most teeth, but it did!"

"Capsules of composite come in a handy box with a magnetic closure."

Evaluation Summary:

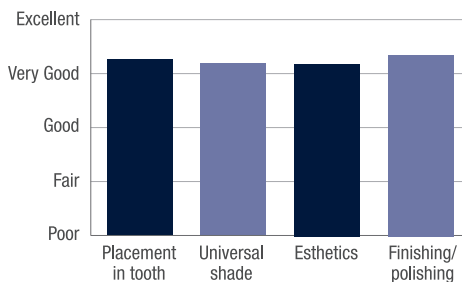

Compared to Competitive Products:

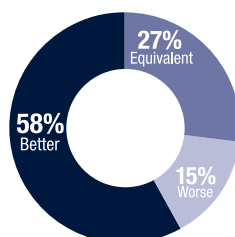

Percentage of Consultants Who Would:

